



Making **cutting-edge**  
medical equipment  
**affordable**

BEDS &  
ACCESSORIES



OVERBED  
TABLES



EVACUATION



BED RAIL  
BUMPERS



WASHABLES



HOISTS &  
TRANSFERS



TOILETING




WHEELCHAIRS



WALKERS &  
ROLLATORS




FALL  
PREVENTION

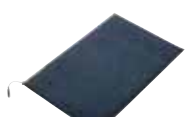


HIP  
PROTECTION




FOAM PRESSURE  
CARE

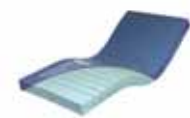

CUSHION PRESSURE  
CARE

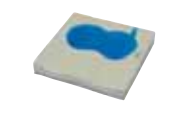

ALTERNATING  
PRESSURE CARE

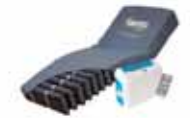

HYBRID PRESSURE  
CARE



SUCTION &  
NEBULISATION





Alerta Medical are a world leader in the design, manufacture and supply of innovative and trusted durable medical equipment, delivering solutions to proactively maximise user safety and comfort.

**Distributor Commitment**

Alerta Medical supply exclusively through a network of distributor partners across the UK and in over 25 countries worldwide. The Alerta brand has become renowned for innovation, quality and reliability, which combined with our industry leading distributor support, provides a unique advantage to our distributors through enhanced brand image and increased profitability.

**Delivering Frontline Solutions**

Alerta Medical understand care providers need dependable and robust equipment to fulfil the rigorous requirements of care and meet high expectations. Alerta Medical was established through a belief that prevention is better than cure and everyone in care should enjoy the benefits of a proactively safer care environment.

**Alerta Medical Mission**

Proactively improving the comfort and quality of life for people in care by making cutting-edge medical equipment affordable to all through strategic distributor partnerships. Everyone in care will be kept comfortable, safe and secure when protected by Alerta.

**Quality Accreditation**

Our medical equipment is manufactured in ISO13485 approved production facilities with thorough quality control procedures, providing distributors with total assurance, strongly backed by our product expertise, delivery capability and customer support, underpinned with a clear vision for the future.

**Alerta Medical Vision**

To be the partner of choice for companies operating in the healthcare market, by building transparent relationships, unrelenting innovation and the pursuit of excellence.

**Proven Experience**

An independent and family owned company, Alerta Medical builds on a legacy of three generations, dedicated to improving clinical outcomes and user experience at the frontline of healthcare.

**We look forward to serving you!**



Certificate Number:  
13617-QMS-002  
ISO9001



**Head Office**

Alerta Medical Ltd  
4 Symington Place  
Riverside Business Park  
Irvine  
KA11 5DE  
United Kingdom



### UNBEATABLE QUALITY

Clinically proven over many years, our products are designed to consistently perform in the healthcare environment – **enhancing your brand image as a distributor.**



### MANUFACTURER PRICES

As the first step in the supply chain, we offer an attractive pricing structure – **increasing your profitability as a distributor.**



### PRODUCT EXPERTISE

As a **global designer and manufacturer** we provide in depth product knowledge whenever you need it.



### PERSONAL SERVICE

As an approved distributor partner you have full access to our **industry leading product range and personal support** from our team.



### LARGE STOCKHOLDING

Our large distribution centre enables us to hold **vast stocks for next day delivery** to you, or direct to your customer.





# Contents

## Beds & Accessories

Lomond Standard Bed	6
Lomond Low Bed	6
Lomond Floor2 Bed	7
Lomond Bariatric Bed	7
Lomond Rotate-Stand-Up Bed	8
Overbed Tables	9
Bed Accessories	10
Crash Mat & Evacuation	11
Bed Safety Bumpers	12-13

## Washables

Bed Pads	14
Chair Pads	14
Linen Bags	15
Bibs	15

## Toileting, Moving & Handling

Powerlifter Hoist Mini 135	18
Powerlifter Hoist Midi 150	18
Powerlifter Hoist Maxi 175	19
Powerlifter Hoist Maxi 175-e	19
MoveAssist 200 Transfer Aid	20
Powerlifter Stand-Aid 200	20
Aqua Chair	21
Wheelchairs	22
Shower Chairs & Commodes	23
Toileting Aids	24
Walkers & Rollators	25
Slide Sheets	26

## Foam, Gel & Static Systems

Product Overview	28
Sensaflex Underlay	29
Foam Mattress Sensaflex 500	29
Foam Mattress Sensaflex 1000	30
Foam Mattress Sensaflex 3000	30
Foam Mattress Sensaflex 4000	31
Foam Mattress Sensaflex 4500	32
Foam Mattress Sensaflex 5000	32
Foam Cushions	33
Gel Cushion - GelCube	34
Gel Cushions	35

5-13	Static Air - Cushion	35
	Static Air - Bootees	35

## Fall Prevention

37-49	Product Overview	38
	Deluxe Alertamat	39
	Wood Alertamat	40
	Alertamat Crash Mat	40
	Floor+ Alertamat	41
	Bed Alertamat	42
	Chair Alertamat	43
14-15	Detect Motion PIR Sensor	44
	Alerta Nurse Call Buttons	45
	WetSense Mats	46
	Accessories & Plug Identification	47
	Safehip Hip Protection	48-49

## 17-26

## Alternating Air Systems

51-75	Product Overview	52
	Air Mattress Bubble 2 Overlay	53
	Air Mattress Pearl Overlay	54-55
	Air Mattress Emerald2 Overlay	56-57
	Air Mattress Emerald Auto Overlay	58-59
	Air Mattress Sapphire2 Replacement	60-61
	Air Mattress Ruby2 Replacement	62-63
	Air Mattress Ruby Auto Replacement	64-65
	Air Mattress Bariatric2 Replacement	66-67
	Sensaflo Hybrid Mattress	68-69
	Sensaflo Hybrid Cushion	70-71
	Partner Air Cushion	72
	Mobile Battery Air Cushion	73
27-35	Product Selector Guide	74
	Bespoke Branding	75

## Suction Machines & Nebulisers

76-77	Suction Machines, Mains & Battery	76
	Nebulisers	77

## Diagnostics

78	Thermometers, Non-Contact	78
	Thermometers, Contact	78
	Blood Pressure Monitor	78
	Oximeter	78





# Beds & Accessories

# Care Beds



ALSO AVAILABLE IN WALNUT



20CM EXTENSION AVAILABLE

The **Alerta Lomond Standard Bed** provides the perfect balance between budget, functionality and quality, and comes with the added benefit of Trendelenburg and reverse Trendelenburg tilting feature.

**ALT-LB-SW - Lomond Standard Bed, Walnut**


**ALT-LB-SO - Lomond Standard Bed, Oak**

## KEY FEATURES AND BENEFITS

- > 10 function lockable handset
- > Trendelenburg and reverse Trendelenburg tilting feature
- > Knee brake adjustment
- > Click stop lower leg raising elevation
- > Electrically adjustable backrest
- > Four robust braking 5" castors
- > Four-section mattress support with mattress grips
- > Extra high side rails available
- > 20cm bed length extension kit available
- > Tool-free assembling and dismantling
- > Transport system for storage and transporting

## TECHNICAL SPECIFICATION

- > Max Safe Working Load: 200kg
- > Height Adjustment: 40 - 80cm
- > Outer Dimensions: 102 x 213cm

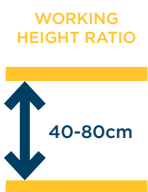

ALSO AVAILABLE IN OAK



20CM EXTENSION AVAILABLE

The **Alerta Lomond Low Bed** provides the ultimate balance between user care, functionality and quality. It has all the features of the Lomond Standard and comes with the added benefit of Trendelenburg and reverse Trendelenburg tilting feature. While being a very low bed, the Lomond Low also has the unique ability of its range extending from 20cm to 80cm high.

**ALT-LB-LW - Lomond Low Bed, Walnut**


**ALT-LB-LO - Lomond Low Bed, Oak**

## KEY FEATURES AND BENEFITS

- > 10 function lockable handset
- > Trendelenburg and reverse Trendelenburg tilting feature
- > Knee brake adjustment
- > Click stop lower leg raising elevation
- > Electrically adjustable backrest
- > Four robust braking 3" castors
- > Four-section mattress support with mattress grips
- > Extra high side rails available
- > 20cm bed length extension kit available
- > Tool-free assembling and dismantling
- > Transport system for storage and transporting

## TECHNICAL SPECIFICATION


- > Max Safe Working Load: 250kg
- > Height Adjustment: 20 - 80cm
- > Outer Dimensions: 102 x 213cm

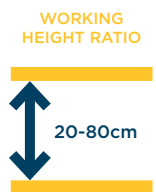

All Lomond Beds have a 3 year warranty on the bed frame and 2 years on the motors and electrics. The beds can be dismantled for storage or transporting and are available in Oak and Walnut wood finish.

# alerta<sup>®</sup> lomond floor2

ALSO AVAILABLE IN WALNUT

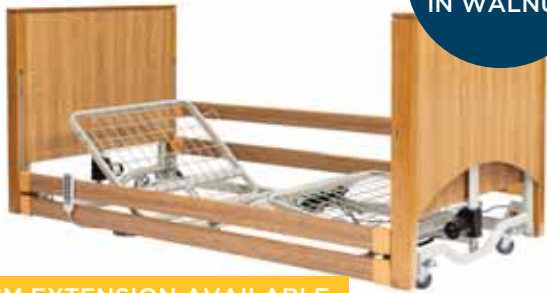

20CM EXTENSION AVAILABLE

The **Alerta Lomond Floor2 Bed** is a unique ultra-low bed for pioneering user safety and care. It has all the features of the Lomond Low and comes with the added benefit of its range going down to just 9cm from the ground, **and auto regression backrest**. This facilitates the most specific and complex care requirements, while also reducing impact if the user rolls out of bed.


- ALT-LB-FW - Lomond Floor Bed, Walnut
- ALT-LB-FO - Lomond Floor Bed, Oak

## KEY FEATURES AND BENEFITS

- > 10 function lockable handset
- > Trendelenburg and reverse Trendelenburg tilting feature
- > Auto regression backrest
- > Knee brake adjustment
- > Click stop lower leg raising elevation
- > Electrically adjustable backrest
- > Four robust braking 3" castors
- > Four-section mattress support with mattress grips
- > Extra high side rails available
- > 20cm bed length extension kit available
- > Tool-free assembling and dismantling
- > Transport system for storage and transporting


## TECHNICAL SPECIFICATION

- > Max Safe Working Load: 200kg
- > Height Adjustment: 9 - 70cm
- > Outer Dimensions: 102 x 213cm

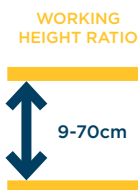

# alerta<sup>®</sup> lomond bari

ALSO AVAILABLE IN OAK

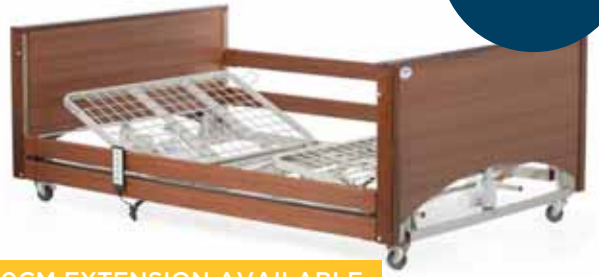

20CM EXTENSION AVAILABLE

The **Alerta Lomond Bariatric 120cm Wide Bed** provides the ultimate balance between user care, functionality and quality. It has all the features of the Lomond Standard and has the added ability of its range extending from 20cm to 80cm high, and a 120cm wide platform.

- ALT-LB-BW - Lomond Bariatric Bed, Walnut
- ALT-LB-BO - Lomond Bariatric Bed, Oak

## KEY FEATURES AND BENEFITS

- > 10 function lockable handset
- > Trendelenburg and reverse Trendelenburg tilting feature
- > Knee brake adjustment
- > Click stop lower leg raising elevation
- > Electrically adjustable backrest
- > Four robust braking 3" castors
- > Four-section mattress support with mattress grips
- > Extra high side rails available
- > 20cm bed length extension kit available
- > Tool-free assembling and dismantling
- > Transport system for storage and transporting

## TECHNICAL SPECIFICATION

- > Max Safe Working Load: 250kg
- > Height Adjustment: 20 - 80cm
- > Outer Dimensions: 132 x 213cm





NEW



The **Alerta Lomond Rotate-Stand-Up Bed** is a powered rotation bed to support users with unassisted transfers between lying in bed and standing, providing enhanced freedom, safety and care for the user.

The user can independently control the bed with the handset to transition from lying in bed to the seated position, and then rotating and standing up from the bed. The hand rails provide additional safety and assurance for the user when rotating and transferring in or out of the bed.

It has a 3 year warranty on the bed frame and 2 years on the motors and electrics. The bed can be dismantled for storage or transporting and is available in Oak wood finish.

The Alerta Lomond Rotate bed also features a battery back-up ensuring the bed never gets stuck in a sitting position in the event of a power cut.

## PRODUCT CODES

**ALT-LB-RO - Lomond Rotate-Stand-Up Bed, Oak**


## KEY FEATURES AND BENEFITS

- > 12 function lockable handset
- > Knee brake adjustment
- > Electrically adjustable backrest
- > Four robust braking 4" castors
- > Four-section mattress support with mattress grips
- > Tool-free assembling and dismantling
- > Memory foam articulated mattress included
- > Hand rails to assist onto and off the bed in chair position
- > One button function from lying to sitting
- > Battery back-up in case of mains failure

## TECHNICAL SPECIFICATION

- > Max Safe Working Load: 250kg
- > Height Adjustment: 40 - 80cm
- > Outer Dimensions: 102 x 213cm
- > Mattress Width: 82cm

WORKING  
HEIGHT RATIO

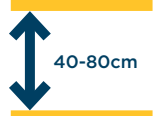



# Overbed Tables



Oak without castors

Walnut with castors

The **Alerta Overbed Table** provides easy access to everyday conveniences for users seated in a bed, chair or wheelchair.

It is finished in oak or walnut and comes with or without castors. The table has a laminate top with two adjusting points which allows both height and tilt adjustment and comes with raised edges to prevent items slipping off.

## TECHNICAL SPECIFICATION

- > Max weight capacity: 10kg
- > Height adjustment range: 64-97cm
- > Table top dimensions: 61 x 40cm
- > Max Tilt angle: 34°

## PRODUCT CODES

ALT-7510 - Alerta Oak Overbed Table, No Castors  
 ALT-751W - Alerta Walnut Overbed Table, No Castors

ALT-7510C - Alerta Oak Overbed Table, With Castors  
 ALT-751WC - Alerta Walnut Overbed Table, With Castors

## KEY FEATURES AND BENEFITS



Castor option

Height & tilt adjustment

Won't topple over even with 10kg on top!



Oak

Walnut

The **Alerta Gas Lift Overbed Table** features an easy lift handle to adjust the table height with minimal effort. It is finished in oak or walnut and comes with brakes on the two rear castors. The table has a laminate top and comes with either a straight leg base or a wheelchair accessible base design.

## TECHNICAL SPECIFICATION

- > Straight Leg or Wheelchair access base
- > Braking castors
- > Max weight capacity: 15kg
- > Height adjustment range: 80-105cm
- > Table top dimensions: 80 x 40cm

## PRODUCT CODES

ALT-GLS-O - Alerta Oak Gas Lift Overbed Table, Straight Leg Base  
 ALT-GLS-W - Alerta Walnut Gas Lift Overbed Table, Straight Leg Base  
 ALT-GLW-O - Alerta Oak Gas Lift Overbed Table, Wheelchair Leg Base  
 ALT-GLW-W - Alerta Walnut Gas Lift Overbed Table, Wheelchair Leg Base

## KEY FEATURES AND BENEFITS



Rear castors

Easy lift height adjustment

Won't topple over even with 15kg on top!


# Bed Accessories

ALT-LB-EK90

LOMOND BED 20CM  
EXTENSION KIT, 90CM WIDE

ALT-LB-EK120

LOMOND BED 20CM  
EXTENSION KIT, 120CM WIDE

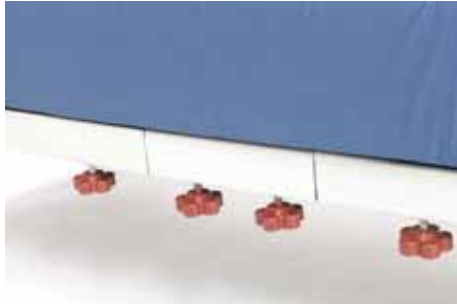

ALT-66

ALERTA FOAM MATTRESS 20CM  
EXTENSION INSERT, 90 X 14CM

ALT-66B

ALERTA FOAM MATTRESS 20CM  
EXTENSION INSERT, 120 X 16CM



ALT-LB-LSRW

LOMOND WALNUT LONG  
SIDE RAILS, SET OF 4

ALT-LB-LSRO

LOMOND OAK LONG  
SIDE RAILS, SET OF 4



ALT-LB-HSRW

LOMOND WALNUT  
SINGLE HIGH SIDE RAIL

ALT-LB-HSRO

LOMOND OAK  
SINGLE HIGH SIDE RAIL



ALT-LB-WP

LOMOND BED  
WALL PROTECTOR



ALT-LB-AL

LOMOND BED  
ASSISTING LEVER



ALT-LB-LP

LOMOND BED LIFTING POLE  
> Comes with strap and handle



ALT-LB-MT

LOMOND BED  
MOVING TROLLEY



NEW


# Crash Mat & Evacuation



## ALT-103 - Folding Crash Mat

The **Alerta Crash Mat** is designed to be placed next to the bed to minimise injury if the user should fall from bed. The mat incorporates triple fold design for easy storage and carry handles for ease of transport.

Alerta Crash Mat is quick and easy to position and benefits from a PU wipe clean surface. When unfolded the mat is 190cm long x 60cm wide x 6cm high.

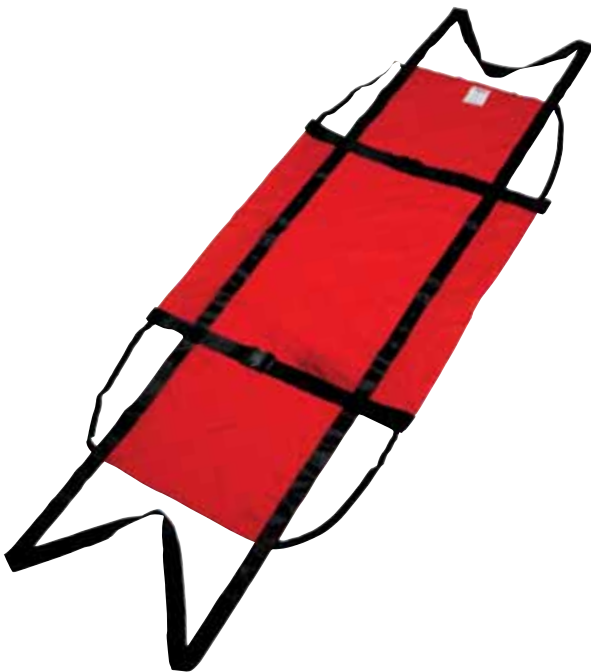

## ALT-E001 - Evacuation Sheet

The **Alerta Evacuation Sheet** fixes permanently to the underside of a flexible foam mattress by means of strong elastics. It has safety belts to hold the user on the mattress and pulling straps at each end.

The combination of slippery sheet and flexible mattress will allow controlled and easy access through doorways and down stairs.

Max user weight: 350kg



Wall mounted zip storage bag



## ALT-E002 - Evacuation Sledge

The **Alerta Evacuation Sledge** is used for transporting a user to safety during an evacuation. The user is placed on the sledge, secured with the chest and leg straps and pulled to safety using the straps at each end.

The flexible sledge and slippery underneath will allow controlled and easy access through doorways and down stairs.

Max user weight: 350kg

# Bed Safety Bumpers

**Alerta Medical** bed protection is made with PU (Polyurethane Fabric) instead of PVC due to risks of Phthalates with PVC. PU is safer and higher quality, stretchable, breathable and serves a longer lifetime.



## ALT-112 - 2 Bar Rail Bumpers with Netting



The Alerta Bumpers with Net Inserts are designed to zip onto side rails, providing protection for users in bed.

Net inserts are centred between the bed rails to improve visibility and comfort for the user.

Suitable for profiling beds with 2 side rails.

200cm strong long-life stretch PU fabric.

**Also available in 20cm extended - code: ALT-115  
220cm strong long-life stretch PU fabric.**

## ALT-113 - 3 Bar Rail Bumpers with Netting



The Alerta Bumpers with Net Inserts are designed to zip onto side rails, providing protection for users in bed.

Net inserts are centred between the bed rails to improve visibility and comfort for the user.

Suitable for profiling beds with 2 side rails and an additional 3rd high side rail attached.

200cm strong long-life stretch PU fabric.

## ALT-114 - Full Length Bed Rail Bumpers



The Alerta Full Length Bed Rail Bumpers are designed to fit onto standard height wooden side rails, providing protection for users in bed.

Suitable for profiling beds with 2 side rails.

200cm strong long-life stretch PU fabric.

**BREATHABLE | IMPERMEABLE | STRETCH | ANTIBACTERIAL | FR**



NEW



### ALT-118 - 2 Bar Rail Anti-Entrapment Bumpers with Netting

Alerta Anti-Entrapment Bumpers benefit from being connected with a mesh base which sits on the mattress platform of the bed underneath the mattress, forming a safety net to help prevent limb entrapment.

The Alerta Anti-Entrapment Bumpers with Net Inserts are designed to velcro onto side rails, providing protection for users in bed. Net inserts are centred between the bed rails to improve visibility and comfort for the user.

Suitable for profiling beds with 2 side rails.



**ALT-116 - Alerta Safety Side Wedge Set** is designed to protect those in care and help keep them comfortably positioned in bed. Made from high density foam with a wipe clean PU cover, it has a draw sheet and two built-in straps that secure the wedge set to the bed. 2 sets can be used down the length of the bed for extra protection if required.

#### KEY FEATURES & BENEFITS

- > Safety wedge
- > Draw sheet included
- > Wipe clean fabric
- > Straps for securing to the bed
- > 99cm wide x 94cm long x 20cm high

#### CARE INSTRUCTIONS

- > Wipe down with a mild detergent or sterilising fluid
- > **Do not use bleach or biological cleaners**
- > **Do not iron**
- > **Do not use fabric softeners**



**ALT-117 - Alerta Head & Foot Board Protector** is designed to protect those in care from hard head and foot boards on their care bed. Made from a wipe clean fabric, each bumper is tapered to allow for tracking side rails and is easy to fit. Use in conjunction with side bumpers to create a comprehensive bed safety system.

#### KEY FEATURES & BENEFITS

- > Safety padding
- > 105cm wide
- > 58cm high
- > Wipe clean fabric
- > Tapered for neat side rail use
- > Use with Alerta side bumpers for all round protection

#### CARE INSTRUCTIONS

- > Wipe down with a mild detergent or sterilising fluid
- > **Do not use bleach or biological cleaners**
- > **Do not iron**
- > **Do not use fabric softeners**

**BREATHABLE | IMPERMEABLE | STRETCH | ANTIBACTERIAL | FR**

# Washables

## ALERTA WASHABLE BED PAD

- > 300 wash guarantee
- > Tough durable construction
- > Fully waterproof (integral barrier)
- > Generously sized tuck-in flaps
- > Machine washable
- > Machine dryable (low heat only)

DBP/2

2.1 LITRE ABSORBENCY, 72X90CM

DBP/3

3.2 LITRE ABSORBENCY, 86X90CM

DBP/D

4 LITRE ABSORBENCY, 86X140CM, DOUBLE



Packaged for retail

## ALERTA WASHABLE CHAIR PAD

- > 300 wash guarantee
- > Tough, durable construction
- > Fully waterproof (integral barrier)
- > Generously sized (50cm x 60cm)
- > Machine washable
- > Machine dryable (low heat only)
- > 1.0 litre absorbency capacity
- > Available in 5 colours



Packaged for retail

DSP/BR

BROWN

DSP/GR

DARK GREEN

DSP/LC

WINE

DSP/PB

PETROL BLUE

DSP/TA

TARTAN

### ALERTA LINEN BAG

- > Tough polyester construction
- > Heavy duty draw cord with toggle
- > Machine washable at 60-70°C
- > Machine dryable (low heat only)
- > Available in 5 colours

- L/BAG/W WHITE
- L/BAG/B BLUE
- L/BAG/Y YELLOW
- L/BAG/G GREEN
- L/BAG/R RED



### ALERTA BIB

- > 300 wash guarantee
  - > Tough durable construction
  - > Fully waterproof (integral barrier)
  - > Machine washable
  - > Machine dryable (low heat only)
  - > Generous 45cm width\*
  - > Standard - 70cm
  - > Long - 90cm
  - > 2 stud fastening positions at neck
- \* all sizes are approximate



Packaged for retail




- ALT-BIB-PS ALERTA ADULT BIB - TARTAN PINK - 70CM
- ALT-BIB-BL ALERTA ADULT BIB WITH CRUMB TRAY - TARTAN BLUE - 90CM
- ALT-BIB-PB ALERTA ADULT BIB WITH CRUMB TRAY - PLAIN BLUE - 90CM
- ALT-BIB-PG ALERTA ADULT BIB WITH CRUMB TRAY - PLAIN GREEN - 90CM
- ALT-BIB-PW ALERTA ADULT BIB WITH CRUMB TRAY - PLAIN WINE - 90CM
- ALT-BIB-MA ALERTA ADULT BIB WITH CRUMB TRAY - MALE DESIGN - 90CM
- ALT-BIB-FM ALERTA ADULT BIB WITH CRUMB TRAY - FEMALE DESIGN - 90CM



ALT-BIB-PS



ALT-BIB-BL

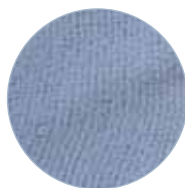

ALT-BIB-PB



ALT-BIB-PG




ALT-BIB-PW



ALT-BIB-MA



ALT-BIB-FM

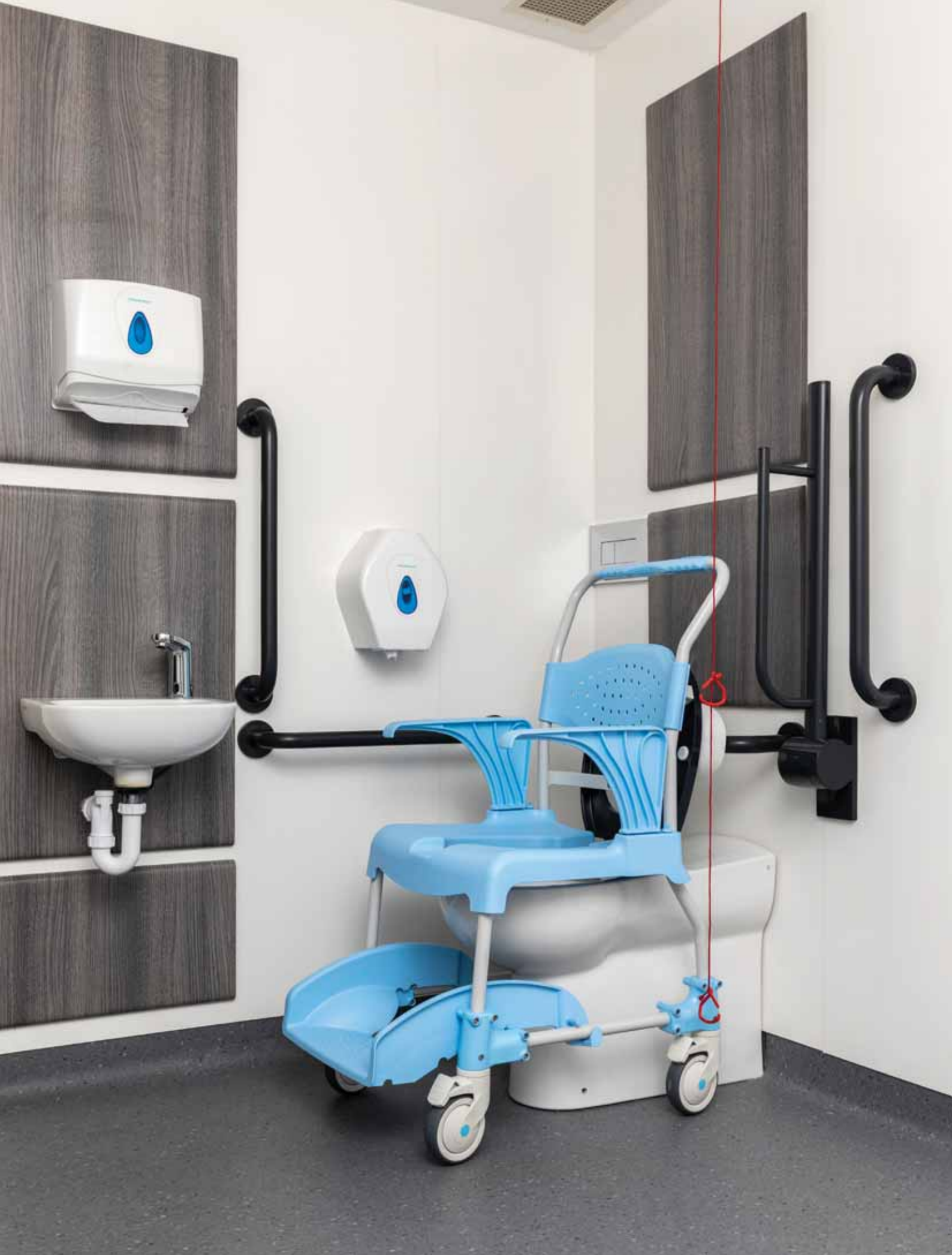





Toileting,  
Moving & Handling

# Patient Transfer Aids



-  Safe working load  
**135kg**
-  Maximum lift height  
**1350mm**
-  Minimum lift height  
**330mm**
-  Overall width  
**530mm**
-  Overall length  
**1130mm**

## ALERTA POWERLIFTER MINI 135

**ALT-PL135 - Alerta Powerlifter Mini 135** – a transfer hoist for handling and transferring users. Suitable for users up to 135 kg. **Alerta Powerlifter** hoists are designed to be used in hospital, nursing and care home environments.

They provide security, functionality and outstanding reliability, easy to operate, easy to move and incorporate anti-entrapment techniques. **Alerta Powerlifter** hoists also give users a very restful and calming lift experience with a thermoplastic rubber coated, anti-swing spreader bar.

The wide range of features make the **Alerta Powerlifter** hoists a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

### TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Spreader bar with 360° rotation and anti-swing
- > Braked rear castors
- > Disassembles into two parts



-  Safe working load  
**150kg**
-  Maximum lift height  
**1470mm**
-  Minimum lift height  
**340mm**
-  Overall width  
**540mm**
-  Overall length  
**1230mm**

## ALERTA POWERLIFTER MIDI 150

**ALT-PL150 - Alerta Powerlifter Midi 150** – a transfer hoist for handling and transferring users. Suitable for users up to 150 kg. **Alerta Powerlifter** hoists are designed to be used in hospital, nursing and care home environments.

They provide security, functionality and outstanding reliability, easy to operate, easy to move and incorporate anti-entrapment techniques. **Alerta Powerlifter** hoists also give users a very restful and calming lift experience with a thermoplastic rubber coated, anti-swing spreader bar.

The wide range of features make the **Alerta Powerlifter** hoists a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

### TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Spreader bar with 360° rotation and anti-swing
- > Braked rear castors
- > Disassembles into two parts



-  Safe working load  
**175kg**
-  Maximum lift height  
**1590mm**
-  Minimum lift height  
**450mm**
-  Overall width  
**630mm**
-  Overall length  
**1360mm**

## ALERTA POWERLIFTER MAXI 175

**ALT-PL175 - Alerta Powerlifter Maxi 175** – a transfer hoist for handling and transferring users. Suitable for users up to 175 kg. **Alerta Powerlifter** hoists are designed to be used in hospital, nursing and care home environments.

They provide security, functionality and outstanding reliability, easy to operate, easy to move and incorporate anti-entrapment techniques. **Alerta Powerlifter** hoists also give users a very restful and calming lift experience with a thermoplastic rubber coated, anti-swing spreader bar.

The wide range of features make the **Alerta Powerlifter** hoists a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

### TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Spreader bar with 360° rotation and anti-swing
- > Braked rear castors
- > Disassembles into two parts



-  Safe working load  
**175kg**
-  Maximum lift height  
**1590mm**
-  Minimum lift height  
**450mm**
-  Overall width  
**630mm**
-  Overall length  
**1360mm**

## ALERTA POWERLIFTER MAXI 175-E

**The Alerta Powerlifter Maxi 175-E** is a transfer hoist with electric legspread, used for handling and transferring users. Suitable for users up to 175kg. **Alerta Powerlifter** hoists are designed to be used in hospital, nursing and care home environments.

They provide security, functionality and outstanding reliability, easy to operate, easy to move and incorporate anti-entrapment techniques. **Alerta Powerlifter** hoists also give users a very restful and calming lift experience with a thermoplastic rubber coated, anti-swing spreader bar.

The wide range of features make the **Alerta Powerlifter** hoists a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

### TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Spreader bar with 360° rotation and anti-swing
- > Braked rear castors
- > Disassembles into two parts



## ALERTA MOVEASSIST 200

**ALT-ST200 - Alerta MoveAssist** - a transfer assist unit which makes transfer functions of all types quick and easy, requiring minimal caregiver assistance. **Alerta MoveAssist** promotes user engagement, enhancing mobility.

Each unit is equipped with a crossbar where users can grasp and pull their self up into a standing position using their own strength. The padded split seats swing out allowing the user to stand up and board the unit. The caregiver then flips the split seats back into place to form a secure and comfortable seat for the user to sit on for the remainder of the transfer.

A user who qualifies to use the **Alerta MoveAssist** must have enough leg and lower body strength to stand up and remain in the standing/ sitting position. Adequate arm strength is required if the user must use the crossbar. Users who satisfy these criteria and have difficulty walking or present a potential for falling down will find the **Alerta MoveAssist** a useful and safe transfer device. For users who lack these requirements, a standaid hoist is recommended.

### TECHNICAL SPECIFICATION

- > Adjustable legs for ease of manoeuvring
- > Ideal for bed-chair-wheelchair-commode transfers
- > Perfect for room to room transfers
- > MoveAssist can be wheeled directly up to a toilet, offering enhanced dignity
- > Braked rear castors
- > Max user weight 200Kgs
- > Sling hooks for extra security
- > Leg height clearance: 10cm



Safe working load  
**200kg**



Maximum lift height  
**1580mm**



Minimum lift height  
**1140mm**



Overall width  
**540 to 990mm**  
(extending legs)



Overall length  
**1130mm**

## ALERTA STAND 200 POWERLIFTER

**ALT-PLS200 - Alerta Stand 200 Powerlifter** - a lifting hoist for assisting users with standing from a seated position. Suitable for users up to 200kg. Alerta Stand lifter is designed to be used in hospital, nursing and care home environments.

It provides security, functionality and outstanding reliability, easy to operate and easy to move. The Alerta Stand lifter also gives users a very restful and calming lifting experience.

The wide range of features make the Alerta Stand lifter a highly versatile and cost-effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

### TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Braked rear castors
- > Disassembles into two parts
- > Adjustable arm with 5 positions
- > Adjustable legs for ease of manoeuvring
- > Max user weight: 200kg





## SHOWER CHAIR | COMMODE CHAIR | TRANSFER CHAIR | TOILETING CHAIR

### ALT-SC150 ALERTA AQUA SHOWER COMMODE CHAIR

The **Alerta Aqua** chair combines 4 chairs in 1, performing seamlessly as a shower chair, a commode chair and a transfer chair.

While it can be used for showering, transferring and as a mobile commode, the pot can be removed and the chair can also be wheeled to the toilet, making this a very easy to use toileting chair.

The Alerta Aqua is stylish, effortlessly manoeuvrable and manufactured from easy to clean plastic coated steel. A selection of soft padded seat sections are also available as accessories if desired.

### KEY FEATURES AND BENEFITS

- > Five position footrest
- > 4 individually braked wheels
- > Weight capacity of 150kg
- > Removable arms
- > 2 year guarantee



# Wheelchairs



## ALT-1100 – Alerta Car Transit Wheelchair

The **Alerta Car Transit Wheelchair** is crash tested and comes with drop down armrests, swing away detachable footrests, arm pads, lap belt and rear brakes, with 30cm rear wheels and 20cm front wheels, all with puncture proof tyres. Seat is 40cm deep x 43cm wide.

### KEY FEATURES AND BENEFITS

- > Crash tested
- > Drop down table-entry arms
- > Padded armrests
- > Lap belt
- > Swing away detachable footrests
- > Puncture proof tyres
- > Rear wheel brakes
- > Heel loops

### TECHNICAL SPECIFICATION

- > 30cm rear wheels
- > 20cm front wheels
- > Seat size: 40cm deep x 43cm wide
- > Max user weight: 115kg
- > Product weight: 15.5kgs



## ALT-1300 Alerta Self-Propelled Wheelchair

The **Alerta Self-Propelled Wheelchair** comes with self-propelling rear wheels, detachable arms, padded armrests, lap belt, swing away and detachable footrests, heel loops and rear brakes, with 60cm rear wheels and 20cm front wheels, all with puncture proof tyres. Seat is 40cm deep x 45cm wide.

### KEY FEATURES AND BENEFITS

- > Self-propelled rear wheels
- > Wave profile self-propel rim
- > Detachable table-entry arms
- > Padded armrests
- > Lap belt
- > Swing away detachable footrests
- > Puncture proof tyres
- > Rear wheel brakes
- > Heel loops

### TECHNICAL SPECIFICATION

- > 60cm rear self-propelled wheels
- > 20cm front wheels
- > Seat size: 40cm deep x 45cm wide
- > Max user weight: 115kg
- > Product weight: 17kg

# Shower Chairs & Commodes

**ALT-BE002**

## **ALERTA MOBILE SHOWER CHAIR**



This mobile shower chair has 2 rear braked castors. 136kg weight capacity.

Boxed in 2's

**ALT-BE004**

## **ALERTA MOBILE SHOWER CHAIR, ADJUSTABLE HEIGHT**



This mobile shower chair has adjustable height and 2 rear braked castors. 136kg weight capacity.

Seat Height: 44.5-54.5cm

Boxed individually

**ALT-BE005**

## **ALERTA STATIONARY SHOWER CHAIR, ADJUSTABLE HEIGHT**



This stationary shower chair has adjustable height. 136kg weight capacity.

Seat Height: 41.5-53.3cm

Boxed in 4's

**ALT-BE009**

## **ALERTA STACKING COMMODE, ADJUSTABLE HEIGHT**

The Alerta Stacking Commode is manufactured from powder coated tubular steel and comes with a comfortable padded upholstery on the seat and back, which is heat sealed allowing it to be cleaned without the risk of water ingress. The toilet seat is oval shaped with a clip on/clip off design. A rubber buffer sits in the seat to avoid seat abrasion against the frame. The commode holds a 5 litre round potty.

180kg weight capacity.

Seat Height: 45-60cm

Boxed in 3's



**ALT-1200**

## **ALERTA COMMODE & TRANSFER CHAIR**

The Alerta Commode & Transfer Chair comes with drop down armrests, swing away detachable footrests, arm pads, four braked castors and wipeable PVC seat and back. Seat is 53cm deep x 46cm wide.

- > Stable, robust powder coated steel frame
- > Drop down padded armrests
- > Swing away detachable footrests
- > Full length arm pads
- > Push handle for optimal positioning
- > 4 solid steel 5" swivel braked castors
- > Fold-up footrests
- > Wipeable PVC seat and back
- > Max user weight: 110kg



# Toileting Aids

**ALT-BE006**

## **ALERTA RAISED TOILET SEAT, 2 INCHES**

This raised plastic toilet seat increases the toilet seat height by 2 inches for easier toilet access in nursing homes and home care.

200kg weight capacity.

Boxed in 12's



**ALT-BE007**

## **ALERTA RAISED TOILET SEAT, 4 INCHES**

This raised toilet seat increases the toilet seat height by 4 inches for easier toilet access in nursing homes.

200kg weight capacity.

Boxed in 12's



**ALT-BE008**

## **ALERTA RAISED TOILET SEAT, 6 INCHES**

This raised plastic toilet seat increases the toilet seat height by 6 inches for easier toilet access in nursing homes and home care.

200kg weight capacity.

Boxed in 12's



**ALT-BE001**

## **ALERTA TOILET SEAT AID, ADJUSTABLE HEIGHT**

This toilet seat aid is plastic coated and has a large seating area. It features cut-outs at the front and rear for easy cleaning access. The toileting seat simply clicks into position on the frame and rotates when required.

165kg weight capacity.

Seat Height: 48-62cm.

Boxed in 2's



**ALT-BE003**

## **ALERTA PORTABLE TOILET SURROUND, ADJUSTABLE HEIGHT**

Portable toilet seat surround with height extension legs designed to make using the toilet easier for users in nursing homes and home care.

Overall Height: 68-80.7cm

Boxed individually





# Walking Aids



WATCH ROLLATOR  
PRODUCT VIDEO



NEW



NEW



NEW

## ALT-R014

### THREE-WHEEL LIGHTWEIGHT WALKER WITH BASKET, TRAY & BAG

The **Alerta Three-Wheel Lightweight Walker** is equipped with a basket, tray and bag, and features a triangular foldable frame and adjustable handle height. It has robust 8" wheels and a lightweight frame for easy transport and manoeuvrability.

### TECHNICAL SPECIFICATION

- > Handle height: 82-97cm
- > Max user weight: 115kg
- > Product weight: 6.5kg

### KEY FEATURES & BENEFITS

- > Robust lightweight frame
- > Lightweight and foldable for manoeuvrability and transport
- > Handle lock safety feature
- > Spacious basket and storage bag
- > Includes tray with raised edges
- > Robust 8" wheels for indoor and outdoor use

## ALT-R004

### THREE-WHEEL LIGHTWEIGHT WALKER WITH SEAT & BAG

The **Alerta Three-Wheel Lightweight Walker** is equipped with a seat, backrest and storage pockets, and features a triangular foldable frame and adjustable handle height. It has robust 8" wheels and a lightweight frame for easy transport and manoeuvrability.

- > Handle height: 82-99cm
- > Seat height: 50cm
- > Max user weight: 135kg
- > Product weight: 6.25kg

- > Seat with backrest for users to rest
- > Robust lightweight frame
- > Lightweight and foldable for manoeuvrability and transport
- > Handle lock safety feature
- > Storage pockets
- > Robust 8" wheels for indoor and outdoor use

## ALT-R006

### FOUR-WHEEL ALUMINIUM ROLLATOR

The **Alerta Four-Wheel Aluminium Foldable Rollator** is equipped with a seat, bag, and independently adjustable backrest. It boasts a unique patented robust quick-fold lightweight frame for easy transport and manoeuvrability, and also benefits from quick release adjustable height handles.

- > Handle height: 79-96cm
- > Seat height: 48cm
- > Max user weight: 136kg
- > Product weight: 7.35kg

- > Seat with independently adjustable backrest for users to rest
- > Robust lightweight aluminium frame
- > Unique patented quick-fold design for easy transport and manoeuvrability
- > Handle lock safety feature
- > Spacious storage bag
- > Robust wheels for indoor and outdoor use

# Slide Sheets



**Packaged for retail**

**ALT-SS - Alerta Slide Sheets** boast superior slide performance properties and are designed to facilitate nurses and carers with safe assisted repositioning and transfer of patients on a support surface.

Alerta Slide Sheets are tubular, reusable, and available in 5 convenient sizes. They are manufactured from a very strong highly durable quality fabric which creates an extremely low friction surface providing superior sliding performance during patient repositioning and transfers whilst protecting against skin shear. Alerta Slide Sheets are not intended for use to lift patients.

Fitted with a brightly coloured hanging loop allowing the sheet to be conveniently stored near the user.

## KEY FEATURES AND BENEFITS

- > Brightly coloured hanging loop
- > Superior slide performance
- > Wash at a maximum of 80°C
- > Do NOT use fabric softener or chlorine
- > Slip hazard; Do not leave on the floor
- > Tubular and reusable

## PRODUCT CODES


<b>ALT-SS-S</b>	<b>Alerta Slide Sheet Small 72 x 70cm</b>
<b>ALT-SS-M</b>	<b>Alerta Slide Sheet Medium 122 x 70cm</b>
<b>ALT-SS-L</b>	<b>Alerta Slide Sheet Large 122 x 100cm</b>
<b>ALT-SS-W</b>	<b>Alerta Slide Sheet Wide 145 x 100cm</b>
<b>ALT-SS-XL</b>	<b>Alerta Slide Sheet X-Large 200 x 140cm</b>



Foam, Gel &  
Static Systems

# Product Overview

All foam products are supplied vacuum roll packed

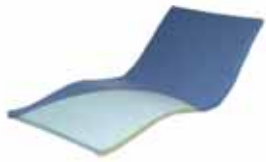

Foam Mattress Underlay



Foam Mattress




Foam Mattress




Memory Foam Mattress



Gel-Topped Memory Foam Mattress

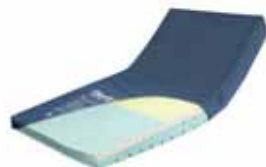

Heel Slope Mattress Memory & Gel Foam

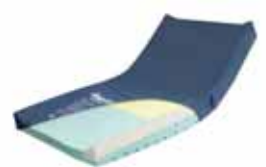

Heel Slope Mattress With Side Supports



Foam Cushion



Memory Foam Cushion





Memory Foam Cushion



Memory Foam Bariatric Cushion





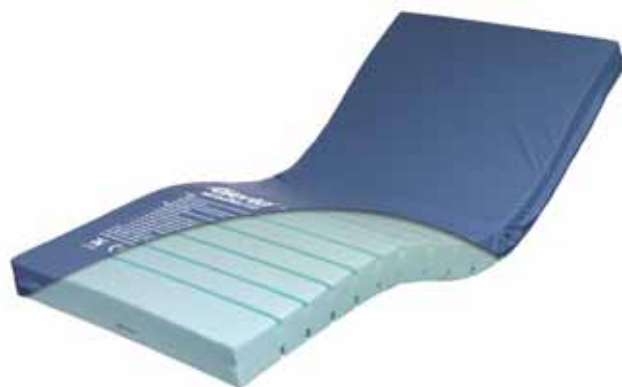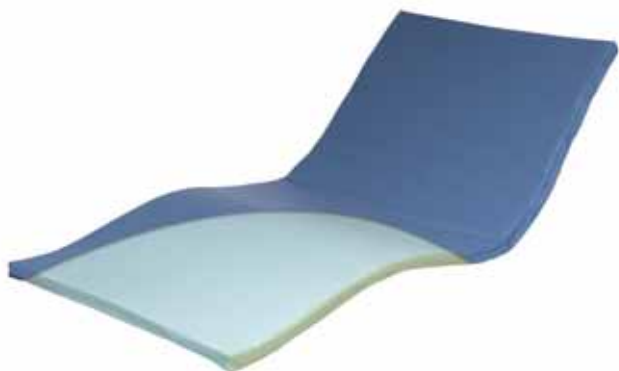

**ALT-UND – Alerta 2 inch Foam Mattress Underlay**  
 – a foam mattress underlay for use underneath an air mattress system in the event of power failure or during the use of CPR. The underlay provides effective comfort and care for users in hospital, nursing and care home environments.

To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

## TECHNICAL SPECIFICATION

### Mattress Specification

- > Mattress dimensions: 198 x 90 x 5 cm (78 x 35.5 x 2")
- > Multi-stretch and vapour permeable PU cover
- > Water resistant and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5

**ALT-500 – Alerta Sensaflex 500** – a medium risk profiling foam mattress built with castellated foam. The sophisticated mattress design provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

## TECHNICAL SPECIFICATION

### Mattress Specification

- > Mattress dimensions: 198 x 90 x 12 cm (78 x 35.5 x 4.8")
- > Max load capacity: 120kg (265lb)
- > Min load capacity: 50kg (110lb)
- > Multi-stretch and vapour permeable PU cover
- > Water resistant and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5



To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

- > Multi-stretch and vapour permeable PU cover
- > Water resistant and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5



**ALT-1000 - Alerta Sensaflex 1000** - a high risk profiling foam mattress built with castellated foam. The sophisticated mattress design provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

## TECHNICAL SPECIFICATION

### Mattress Specification

- > Mattress dimensions: 198 x 90 x 14 cm (78 x 35.5 x 5.5")
- > Max load capacity: 180kg (397lb)
- > Min load capacity: 40kg (88lb)



**ALT-3000 - Alerta Sensaflex 3000** - a high risk profiling memory foam mattress. The sophisticated mattress design provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

## TECHNICAL SPECIFICATION

### ALT-3000 Mattress Specification

- > Mattress dimensions: 198 x 90 x 14 cm (78 x 35.5 x 5.5")
- > Max load capacity: 180kg (397lb)
- > Min load capacity: 40kg (88lb)

### ALT-3000/4 Mattress Specification

- > Mattress dimensions: 198 x 120 x 16 cm (78 x 47 x 6.3")
- > Max load capacity: 200kg (441lb)
- > Min load capacity: 40kg (88lb)



# alerta<sup>®</sup>

## sensaflex4000



LOW  
MINIMUM  
USER  
WEIGHT



**ALT-4000 - Alerta Sensaflex 4000** - a very high risk profiling memory foam mattress built with a high density base, foam centre and a castellated gel memory foam top. The mattress is designed with a section of v-cut memory foam which creates additional air movement for enhanced pressure relieving properties, providing effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

## TECHNICAL SPECIFICATION

### Mattress Specification

- > Mattress dimensions:  
198 x 90 x 16 cm  
(78 x 35.5 x 6.3")
- > Max load capacity: 200kg (441lb)
- > Min load capacity: 30kg (66lb)
- > Multi-stretch and vapour permeable PU cover
- > Water resistant and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5



**All foam mattresses (except 2" underlay) and cushions are shipped vacuum rolled for ease of transport and minimising storage**

# Heel Slope Mattresses




Alerta heel slope mattresses are very high risk profiling memory foam mattresses built with a high density base, memory foam top, and gel-infused memory foam sloping heel section. They are designed with a section of v-cut base foam which creates additional air movement for enhanced pressure relieving properties, providing effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

- > Heel slope: assists pressure redistribution by alleviating pressure from the high risk heel area.
- > SensaGel® heel zone: gel-infused memory foam to prevent pressure sores on the heel.

To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

- > Multi-stretch and vapour permeable PU cover
- > Water resistant and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5

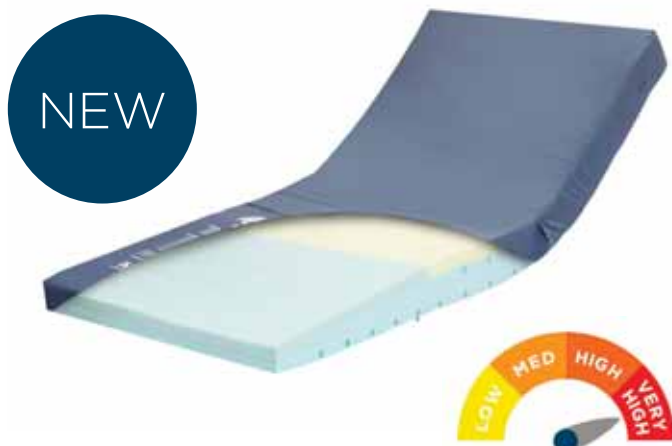

**ALT-4500 - Alerta Sensaflex 4500** - Memory gel-foam heel slope mattress.

**ALT-5000 - Alerta Sensaflex 5000** - Memory gel-foam heel slope mattress with built-in side wedge supports.

- > Side wedge support: prevents the user rolling over the edge of the mattress, reducing risk of falling out of bed. Also provides anti-entrapment by preventing users limbs becoming trapped in side rails.
- > Patent Number: 2209132.6

## TECHNICAL SPECIFICATION

### Mattress Specification

- > Mattress Dimensions: 198 x 90 x 8-14cm (78 x 35.5 x 3.1-5.5")
- > Max Load: 180kg (397lb)
- > Min Load: 40kg (88lb)



## TECHNICAL SPECIFICATION

### Mattress Specification

- > Mattress Dimensions: 198 x 90 x 14-20 cm (78 x 35.5 x 5.5-7.8")
- > Max Load: 180kg (397lb)
- > Min Load: 40kg (88lb)





To launder the cover in a washing machine set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

- > Multi-stretch and vapour permeable PU cover
- > Water resistant and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5

## alerta<sup>®</sup> sensaflex200



**ALT-200 - Alerta Sensaflex 200** – a medium risk foam cushion built with castellated foam. The foam cushion provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

### TECHNICAL SPECIFICATION

- > Cushion dimensions:  
43 x 43 x 10 cm (17 x 17 x 4")
- > Max Load: 150kg (331lb)
- > Min Load: 40kg (88lb)



## alerta<sup>®</sup> sensaflex250



**ALT-250 - Alerta Sensaflex 250** – a high risk memory foam cushion which provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

### TECHNICAL SPECIFICATION

- > Cushion dimensions:  
43 x 43 x 10 cm (17 x 17 x 4")
- > Max Load: 150kg (331lb)
- > Min Load: 40kg (88lb)



## alerta<sup>®</sup> sensaflex300



**ALT-300 - Alerta Sensaflex 300** – a very high risk memory foam cushion built with a high density base, foam centre and a castellated gel memory foam top. The cushion provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

### TECHNICAL SPECIFICATION

- > Cushion dimensions:  
43 x 43 x 10 cm (17 x 17 x 4")
- > Max Load: 150kg (331lb)
- > Min Load: 40kg (88lb)



## alerta<sup>®</sup> sensaflex 350





**ALT-350 - Alerta Sensaflex 350** – a high risk bariatric memory foam cushion which provides effective comfort, care and pressure redistribution for bariatric users in hospital, nursing and care home environments.

### TECHNICAL SPECIFICATION

- > Cushion Dimensions:  
45 x 56 x 12 cm (17.7 x 22 x 5")
- > Max Load: 280kg (617lb)
- > Min Load: 50kg (110lb)




# Gel Cushions

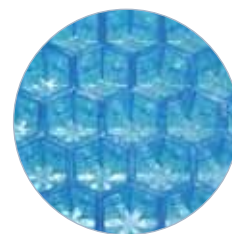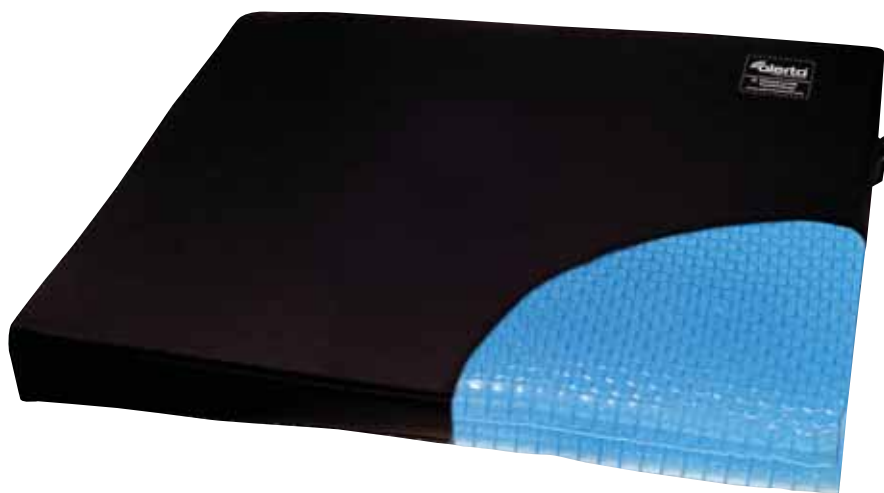

Patented dual honeycomb structure



Intelligent design for balanced pressure distribution

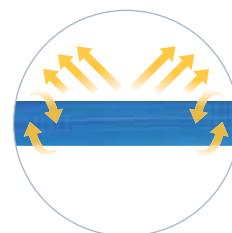

Provides enhanced blood circulation and skin health



Very easy to clean



Packaged for retail




**ALT-204 - Alerta GelCube Cushion** - a patented dual honeycomb air-gel cushion to provide advanced comfort, care and pressure redistribution for users in hospital, nursing and care home environments. It is built with an intelligent dual honeycomb structure to provide optimal blood circulation and skin health, and balanced pressure distribution. The combination of flexible air-gel cells provides a non-slip surface for superior comfort over an extended period of use.

The Honeycomb Gel structure can be washed using soap and warm water up to max 70°C. Do not leave submerged in water for a long period of time and do not squeeze or twist the cushion after washing. Allow to dry naturally on a flat surface. Do not store in direct sunlight.


The bottom of the cover is coated with silicon for non-slip and the top of the cover is made with polyester fabric. The cover can be washed in a washing machine by putting it in a laundry net up to 40°C. Tumble drying is allowed if the cover is in a laundry net.

## KEY FEATURES & BENEFITS


- > Patented dual honeycomb structure
- > Intelligent design for balanced pressure distribution
- > Provides enhanced blood circulation and skin health
- > Cushion dimensions: 46 x 42 x 5 cm (18 x 17 x 2")
- > Waterproof and fire-retardant cover
- > Non-slip polyester cover

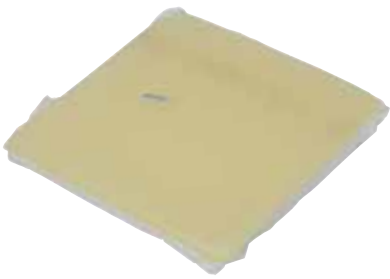

Packaged for retail

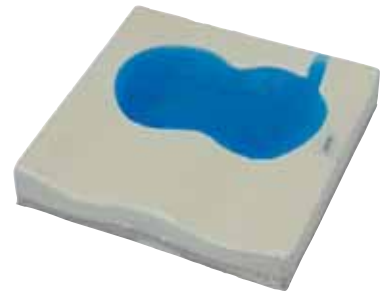

Cover provided

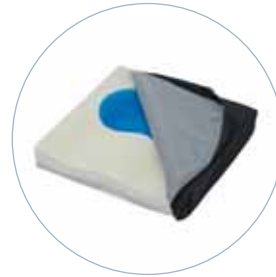

Cover provided



**ALT-201 - Alerta Gel Cushion** - a 43x43x2.5cm gel cushion containing polyurethane SensaGel which acts as a layer of artificial fat, to provide effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments. It can slide and move laterally to effectively reduce shear forces and is guaranteed not to harden over time, neither leak nor dry if accidentally punctured.

The PU cover is fire resistant classified M1, waterproof and vapour permeable. It can be laundered at 90°C. The cushion can be cleaned with a damp cloth. For full cleaning instructions please see user manual.

**ALT-202 - Alerta Gel-Visco Cushion** - a 43x43x8cm gel-visco cushion built with a contoured moulded viscoelastic foam base, connected to a polyurethane SensaGel insert with two densities, to provide effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments. The viscoelastic SensaGel provides maximum relief to bony areas and the shape ensures a good posture while seated and prevents users from sliding forward.

The PU cover is fire resistant classified M1, waterproof and vapour permeable. It can be laundered at 90°C. The cushion can be cleaned with a damp cloth. For full cleaning instructions please see user manual.



**ALT-203 - Alerta Air Cushion with Pump** - an aid in the prevention of pressure ulceration for all users including those assessed as high risk, and for the treatment of existing stage 1 and 2 pressure ulcers. Clinical supervision is advised if skin damage is severe. Once cleaned this product can be reused by another user.



**ALT-210 - Alerta Air Heel Bootees with Strap** - provides effective protection for users by offloading pressure from the heel to minimise the risk of pressure ulceration and injury to the foot, including those assessed as very high risk. If existing pressure related tissue injury is severe then clinical supervision is recommended.







WATCH WIRELESS  
PRODUCT VIDEO

# Fall Prevention

# Product Overview



**alerta**  
deluxe

Premium floor mat,  
Wired and Wireless



**alerta**  
wood

Premium floor  
mat, Wired



Crash Mat,  
Wireless



**alerta**  
floor+

Floor mat, Wired  
and Wireless



**alerta**  
bed

Bed mat system,  
Wired and Wireless



**alerta**  
chair

Chair mat system,  
Wired and Wireless



**alerta**  
detect

PIR motion sensor,  
Wired and Wireless



**alerta**  
button

Nurse call buzzer,  
Wired and Wireless



**alerta**  
wetsensechair

WetSense  
Chair Mat



**alerta**  
wetsensebed

WetSense  
Bed Mat





The **Deluxe Alertamat** is a robust long-life luxury carpet Alertamat; the world's first and only genuine carpeted Alertamat. Available in 5 colours to suit an existing colour scheme, the **Deluxe Alertamat** provides a more refined touch for the care environment.

The **Deluxe Alertamat** is most effective when placed at the bedside or doorway; activating the alarm when a user leaves their safe zone by stepping on the mat. This alerts the caregiver to ascertain the user's safety. Selecting a colour that blends with the existing flooring makes it an ideal solution for situations where a dementia user may otherwise try to avoid the mat.

Both **Wired** and **Wireless** systems can be connected to a nurse call system for centralised activation, or used on their own, independent of a nurse call system to give localised alerts.

AVAILABLE IN 5 COLOURS



- > Genuine non-slip base
- > Luxury carpet finish
- > Choice of 5 colours
- > Long lasting welded seals and strong wiring
- > Large size avoids chance of overstepping
- > Dementia friendly
- > Dimensions: 60 x 90 cm (23.5 x 35.5")



**WIRED DELUXE ALERTAMAT**

**DAM - Wired Deluxe Alertamat** - can be connected to a nurse call system with a 3 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (sold separately).

- > Plug-and-play system
- > System includes:  
1x Deluxe Alertamat and 1x 3 metre cable

The Deluxe Alertamat comes with a unique built-in "Y" cable, providing a spare socket for other Alerta devices to be connected - using only one nurse call point.



Connect to nurse call system or stand-alone alarm monitor

Breakaway cable to reduce trip hazard and prevent damage



**WIRELESS DELUXE ALERTAMAT**

**W-DAM-T - Wireless Deluxe Alertamat** - comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > System includes: 1x Deluxe Alertamat and 1x Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

**Alerta Wireless Transmitter**

- > Battery powered
- > Takes 2x 1.5v AAA batteries



**Alerta Wall Point Receiver**

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries





**WAM-FHOAK - Wood Alertamat** - a robust long-life luxury Alertamat with a wood safety floor finish; the world's first and only genuine wood finish Alertamat. The Wood Alertamat provides a more refined touch for the care environment.

The **Wood Alertamat** is most effective when placed at the bedside or doorway; activating the alarm when a user leaves their safe zone by stepping on the mat. This alerts the caregiver to ascertain the user's safety. The light colour makes it an ideal solution for situations where a dementia user may try to avoid that mat.

It can be connected to a nurse call system with a 3 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (sold separately).

The **Wood Alertamat** comes with a unique built-in "Y" cable, providing a spare socket for other Alerta devices to be connected - using only one nurse call point.

## KEY FEATURES AND BENEFITS

- > Plug-and-play system
- > Genuine non-slip base
- > Luxury wood finish
- > Long lasting welded seals and strong wiring
- > Hygienic, antibacterial easy-clean material
- > Large size avoids chance of overstepping
- > Dementia friendly
- > Dimensions: 55 x 80 cm (21.5 x 31.5")
- > System includes: 1x Wood Alertamat and 1x 3 metre cable



Connect to nurse call system or stand-alone alarm monitor



Breakaway cable to reduce trip hazard and prevent damage

**W-ACM-T - Alertamat Crash Mat with Wireless Nurse Call Alert** - is designed to be placed next to the bed to minimise injury if the user should fall from bed and activates the nurse call system when a fall occurs.

The mat incorporates sensor pads and a Wireless Transmitter within the foam which connects to the Alerta Wall Point Receiver. It can be connected to a nurse call system for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Wall Point Receiver only. The volume can be deactivated for silent in-room operation when connected to a nurse call system.

This helps prevent a user remaining on the mat undetected for an unknown length of time after falling, by alerting the caregiver to provide timely assistance.

The Alertamat Crash Mat is quick and easy to position and benefits from a PU wipe clean surface. When unfolded the mat is 180cm long x 60cm wide x 6cm high. The mat folds in half for easy storage and transportation.

## KEY FEATURES AND BENEFITS

- > Get alerted when a user falls from bed
- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > 95m wireless range from transmitter to receiver
- > Genuine non-slip base
- > Dimensions: 180 x 60cm (71 x 23.5")
- > Thickness: 6cm (2.5")
- > System includes: 1x Alertamat Crash Mat and 1x Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

### Alerta Wireless Transmitter

- > Battery powered
- > Takes 2x 1.5v AAA batteries



### Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries







The **Floor+ Alertamat** is a robust long-life Alertamat and is most effective when placed at the bedside or doorway; activating the alarm when a user leaves their safe zone by stepping on the mat. This alerts the caregiver to ascertain the user's safety. The light colour makes it an ideal solution for situations where a dementia user may try to avoid the mat.

Both **Wired** and **Wireless** systems can be connected to a nurse call system for centralised activation, or used on their own, independent of a nurse call system to give localised alerts.



WATCH WIRELESS  
PRODUCT VIDEO

- > Genuine non-slip base
- > Long lasting welded seals and strong wiring
- > Hygienic, antibacterial easy-clean material
- > Large size avoids chance of overstepping
- > Dementia friendly
- > Dimensions: 60 x 90 cm (23.5 x 35.5")



## WIRED FLOOR+ ALERTAMAT

**PAM - Wired Floor+ Alertamat** - can be connected to a nurse call system with a 3 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (sold separately).

- > Plug-and-play system
- > System includes:  
1x Floor+ Alertamat and 1x 3 metre cable

The Floor+ Alertamat comes with a unique built-in "Y" cable, providing a spare socket for other Alerta devices to be connected - using only one nurse call point.



Connect to nurse call system or stand-alone alarm monitor

Breakaway cable to reduce trip hazard and prevent damage



## WIRELESS FLOOR+ ALERTAMAT

**W-PAM-T - Wireless Floor+ Alertamat** - comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > System includes: 1x Floor+ Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

### Alerta Wireless Transmitter

- > Battery powered
- > Takes 2x 1.5v AAA batteries



### Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries





The **Bed Alertamat** is most effective when concealed under the bed sheet on top of the mattress; activating the alarm when a user leaves their bed. This alerts the caregiver to ascertain the user's safety. The alarm is reset automatically when the user returns to the bed.

Both **Wired** and **Wireless** systems can be connected to a nurse call system for centralised activation, or used on their own, independent of a nurse call system to give localised alerts.

- > Hygienic, antibacterial easy-clean material
- > Dimensions: 25 x 76 cm (10 x 30")



WATCH WIRELESS PRODUCT VIDEO



## WIRED BED ALERTAMAT

**BAM - Wired Bed Alertamat** - can be connected to a nurse call system with a 2 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (supplied as part of the system).

- > Plug-and-play system
- > System includes: 1x Bed Alertamat, 1x Alerta Alarm Monitor and 1x 2 metre cable

### Alerta Alarm Monitor

- > High, low and silent volume control
- > Multiple choice alarm tone
- > Adjustable delay to reduce false alarms
- > Low battery alert
- > Mains or battery powered
- > Takes 1x 9v battery




Battery operated



Mains operated



Alerta Alarm Monitor

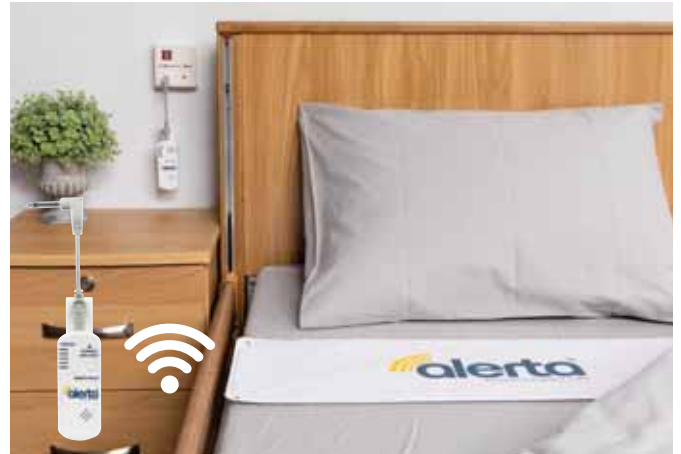

## WIRELESS BED ALERTAMAT

**W-BAM-T - Wireless Bed Alertamat** - comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > 95m wireless range from transmitter to receiver
- > No cables, extra safe and versatile
- > System includes: 1x Bed Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

### Alerta Wireless Transmitter

- > Battery powered
- > Takes 2x 1.5v AAA batteries



### Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries





The **Chair Alertamat** is most effective when concealed under the chair cushion; activating the alarm when a user leaves their chair. This alerts the caregiver to ascertain the user's safety. The alarm is reset automatically when the user returns to the chair.

Both **Wired** and **Wireless** systems can be connected to a nurse call system for centralised activation, or used on their own, independent of a nurse call system to give localised alerts.

- > Hygienic, antibacterial easy-clean material
- > Dimensions: 25 x 38 cm (10 x 15")



WATCH WIRELESS  
PRODUCT VIDEO



## WIRED CHAIR ALERTAMAT

**CAM - Wired Chair Alertamat** - can be connected to a nurse call system with a 2 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (supplied as part of the system).

- > Plug-and-play system
- > System includes: 1x Chair Alertamat, 1x Alerta Alarm Monitor and 1x 2 metre cable

### Alerta Alarm Monitor

- > High, low and silent volume control
- > Multiple choice alarm tone
- > Adjustable delay to reduce false alarms
- > Low battery alert
- > Mains or battery powered
- > Takes 1x 9v battery



Battery operated



Mains operated



Alerta Alarm Monitor



## WIRELESS CHAIR ALERTAMAT

**W-CAM-T - Wireless Chair Alertamat** - comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > 95m wireless range from transmitter to receiver
- > No cables, extra safe and versatile
- > System includes: 1x Chair Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

### Alerta Wireless Transmitter

- > Battery powered
- > Takes 2x 1.5v AAA batteries



### Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries





The **Alerta Detect** uses infrared technology to activate the alarm when body movement is detected within the range of the infrared beam.

The **Alerta Detect** is most effective when mounted at the door or bedside; activating the alarm when a user leaves their safe zone and breaks the infrared beam. This alerts the caregiver to ascertain the user's safety. It is ideal for providing users with the freedom and safety of their own space while still proactively giving the caregiver advanced warning when a user leaves their safe zone.



WATCH WIRELESS PRODUCT VIDEO



## WIRED ALERTA DETECT

**ALT-DET - Wired Alerta Detect** - connects to a nurse call system with a 2 metre cable provided for centralised activation.

### Alerta Detect Motion Sensor

- > Plug-and-play system
- > Mains or battery powered
- > Low battery alert
- > System includes: 1x Alerta Detect motion sensor and 1x 2 metre cable
- > Takes 3x 1.5v AA batteries



## WIRELESS ALERTA DETECT

**ALT-DET - Wireless Alerta Detect** - can be connected to a nurse call system through the Alerta Wall Point Receiver for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Wall Point Receiver only.

- > Wireless connect-and-play system
- > 95m wireless range from transmitter to receiver
- > No cables, extra versatile
- > System includes: 1x Alerta Detect motion sensor (must be used with Alerta Wall Point Receiver - sold separately)

### Alerta Detect Motion Sensor

- > Mains or battery powered
- > Low battery alert
- > Takes 3x 1.5v AA batteries



### Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries





## WIRED ALERTA NURSE CALL BUTTON

**PCL - Wired Alerta Nurse Call Button** - provides the user with reassurance that help is at hand. With a large easy to press button it activates the alarm when a user presses the button. This alerts the caregiver to ascertain the user's need.

It can be plugged directly into the nurse call system or via the spare 'Y' cable that is provided with the wired Deluxe Alertamat, Wood Alertamat and Floor+ Alertamat.

- > 2 metre cable



## ALERTA WALL POINT RECEIVER

**W-WPR - Alerta Wall Point Receiver** - can be connected to a nurse call system for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through this device only. The volume can be deactivated for silent in-room operation when connected to a nurse call system.

The **Alerta Wall Point Receiver** can be wirelessly paired with up to 5 Alerta wireless products and is compatible with the Deluxe Alertamat, Floor+ Alertamat, Bed Alertamat, Chair Alertamat, Alerta Detect and the Alerta Button, Alertmat Crash mat, WetSense Chair Alertamat, and WetSense Bed Alertamat.

It uses radio frequency technology to connect, with a wireless range of up to 95 metres. Once paired with Alerta wireless devices it only transmits to those devices and will not interfere with other Alerta Wall Point Receivers or devices in use nearby.

- > Wireless connect-and-play system
- > No cables, extra versatile
- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries



## WIRELESS ALERTA NURSE CALL BUTTON

**W-NCB - Wireless Alerta Nurse Call Button** - provides the user with reassurance that help is at hand. With a large easy to press button it activates the alarm when a user presses the button. This alerts the caregiver to ascertain the user's need.

It connects to a nurse call system through the **Alerta Wall Point Receiver** for centralised activation and comes with a lanyard, and a holder that can be wall mounted.

- > Wireless connect-and-play system
- > No cables, extra versatile
- > 95m wireless range from transmitter to receiver
- > Takes 1x 12v 23A battery
- > System includes: 1x Alerta Button (must be used with Alerta Wall Point Receiver - sold separately)



The **Alerta Wireless WetSense Bed** and **Chair Alertamats** are designed to inform caregivers when enuresis needs to be addressed. **Alerta WetSense** mats can be used under a sheet on a bed, or on a chair to wirelessly activate the alarm/nursecall system when wetness is detected. This alerts the caregiver to take appropriate action and ensure a comfortable resolution to the individual's needs.

The **Wireless WetSense Alertamats** can be connected to a nursecall system for centralised activation, or used on their own, independent of a nursecall system to give localised alerts.



## WIRELESS WETSENSE CHAIR ALERTAMAT

**W-CAM-WET-T - Wireless WetSense Chair Alertamat**  
- comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > Dimensions: 25 x 38 cm (10 x 15")
- > 95m wireless range from transmitter to receiver
- > System includes: 1x WetSense Chair Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

### Alerta Wireless Transmitter

- > Battery powered
- > Takes 2x 1.5v AAA batteries



### Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries



## WIRELESS WETSENSE BED ALERTAMAT

**W-BAM-WET-T - Wireless WetSense Bed Alertamat**  
- comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > Dimensions: 25 x 76 cm (10 x 30")
- > 95m wireless range from transmitter to receiver
- > System includes: 1x WetSense Bed Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

### Alerta Wireless Transmitter

- > Battery powered
- > Takes 2x 1.5v AAA batteries



### Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries



These products are designed as a monitor alerting the staff to user enuresis. They do not monitor a user leaving the bed or chair.

These mats have a 3 month warranty due to the nature of the substance they are expected to be in regular contact with. The WetSense Alertamats can be replaced on their own and the Alerta Wireless Transmitter and Alerta Wall Point Receiver be reused with a replacement mat.

# Accessories & Plug Identification

## W-FBCT

### ALERTA WIRELESS TRANSMITTER

The Alerta Wireless Transmitter is supplied with the following wireless systems: Deluxe Alertamat, Floor+ Alertamat, Bed Alertamat, Chair Alertamat, Alertmat Crash mat, WetSense Chair Alertamat, and WetSense Bed Alertamat.

It uses radio frequency technology to connect, with a wireless range of up to 95 metres. Once paired with the Alerta Wall Point Receiver it will only transmit to that Wall Point Receiver and will not interfere with other Alerta Wall Point Receivers or devices in use nearby. Requires 2x 1.5v AAA batteries.



## AUAM

### ALERTA ALARM MONITOR

The Alerta Alarm Monitor can be used with all wired Alertamats. It takes 1x 9v battery or can be mains operated.



## MA

### MAINS ADAPTOR

Works with the Alerta Alarm Monitor.



## EXT

### EXTENSION LEAD

Please specify the required length at point of order.



## DPA

### DOUBLE PLUG ADAPTOR

Allows two devices to be connected to one nurse call point.



## Identify your plug type

Most nurse call systems take one of the following plug types:

### SINGLE RING PLUG



### DOUBLE RING PLUG



### WHITE PLUG, 6 GROOVE



### CLEAR PLUG, 6 GROOVE



### CLEAR PLUG, 4 GROOVE



Please ensure that your clients know the plug type they require at point of order to ensure that their Alertamats will function correctly when they receive them.

Plug types can sometimes be difficult to identify. If your client is unsure or if you need assistance, please contact Alerta Medical experts on +44 (0)3452 088 097 or [sales@alertamedical.com](mailto:sales@alertamedical.com)

SCAN HERE FOR  
ALERTA MEDICAL  
PLUG TYPE  
SELECTOR





**HIP PROTECTION**

**CLINICALLY PROVEN**

**SEAMLESS**

**WASHABLE**

### Hip protection helps prevent hip fractures

People who have suffered a fall – or who are at risk of falling – can benefit from effective hip protection. What starts out as an innocent fall may lead to an operation, a long stay in hospital and an increased demand for professional care. In addition, 30% of all people who suffer a hip fracture die within 2 years (1).

A growing number of elderly people indicates a future rise in the occurrence of hip fractures. Worldwide, it is estimated that the number of hip fractures will increase from 1.66 million into 2.6 million in 2025 and 6.26 million in 2050 (2).

Hip protection should be seen as an extra security that contributes to the user's own safety on the same level as seatbelt, lifejacket, or motorcycle helmet. Hip protection reduces the risk of serious injury, saves money and can also improve quality of life by promoting a sense of security, thus allowing the user to remain active, social and reduce the fear of falling.

**SAFEHIP® AirX®** is a 100% textile hip protector with the protector shield based on **AirX®** spacer fabric technology, which is a breathable, skin friendly and soft textile material. The **AirX®** technology is your guarantee for optimal compliance.

**SAFEHIP® AirX®** hip protectors come in three models – a female, male and unisex model in white. The pants are flexible and comfortable to wear, not too warm, and may be used with incontinence pads as a fixation pant. They have been designed to perform like regular underwear, with a fly opening on the male model and lace-look detail on the female. When the gender of the patient is not known in advance, the Unisex brief is suitable for both.

**All SAFEHIP® AirX®** pants are very durable and can be washed in industrial washing machines at up to 95°C/200°F and tumble dried at highest temperatures.

#### Unique features:

- > Patented horse-shoe shaped shield designed to disperse energy from a fall away from the hip joint
- > The shields adapt to the body perfectly, so the hip protector is not noticeable under clothes
- > Excellent breathability (test winner in EMPA test no. 448819)
- > Prevents moisture build-up and feels comfortable against skin
- > Minimises skin irritation such as reddening and itching
- > Pants are very comfortable to wear - no irritating seams
- > High percentage of elastane ensures a perfect fit, which also makes the pants suitable for incontinent users





SAFEHIP®	Cm	Inch	Find your size
S	75-95	30-38	
M	90-110	36-44	
L	100-120	40-48	
XL	110-140	44-56	

Sizes overlap. To find the user's correct size, measure the hips at the widest point. If the hip measurement lies between two sizes, choose the smaller size if the patient has slim thighs, and the larger size if thighs are wide or an incontinence pad will be worn.

**Protector**  
100% AirX® soft textile

Seamless Briefs  
58% Cotton  
36% Polyamide  
6% Elastane



### SAFEHIP AIR-X MALE



Male version with fly

<b>ALT-SAM-S</b>	Air-X Male - S
<b>ALT-SAM-M</b>	Air-X Male - M
<b>ALT-SAM-L</b>	Air-X Male - L
<b>ALT-SAM-XL</b>	Air-X Male - XL

### SAFEHIP AIR-X FEMALE



Female version with flower detail at front

<b>ALT-SAF-S</b>	Air-X Female - S
<b>ALT-SAF-M</b>	Air-X Female - M
<b>ALT-SAF-L</b>	Air-X Female - L
<b>ALT-SAF-XL</b>	Air-X Female - XL

### SAFEHIP AIR-X UNISEX



Neutral design for both sexes

<b>ALT-SAU-S</b>	Air-X Unisex - S
<b>ALT-SAU-M</b>	Air-X Unisex - M
<b>ALT-SAU-L</b>	Air-X Unisex - L
<b>ALT-SAU-XL</b>	Air-X Unisex - XL


### SAFEHIP AIR-X OPEN



Open version

<b>ALT-SAO-S</b>	Air-X Open - S
<b>ALT-SAO-M</b>	Air-X Open - M
<b>ALT-SAO-L</b>	Air-X Open - L
<b>ALT-SAO-XL</b>	Air-X Open - XL

1) IOF - EFFO Jay magaziner, PhD, MSHyg, Eleanor M. Simonsick, PhD, T. Michael Kashner, PhD, J Richard Hebel, PhD, and John E. Konzona, MD Hvidovre Symposium Open Bridge 08/2000  
 2) Cooper et al. Osteoporosis Int. 1992; 2:285-9 & Melton L.J. Bone 1996; 18:121S-5S  
 3) Hindsø K, Lauritzen JB. Intervention study with hip protectors [abstract]. Osteoporosis International, 1998; 8:119.

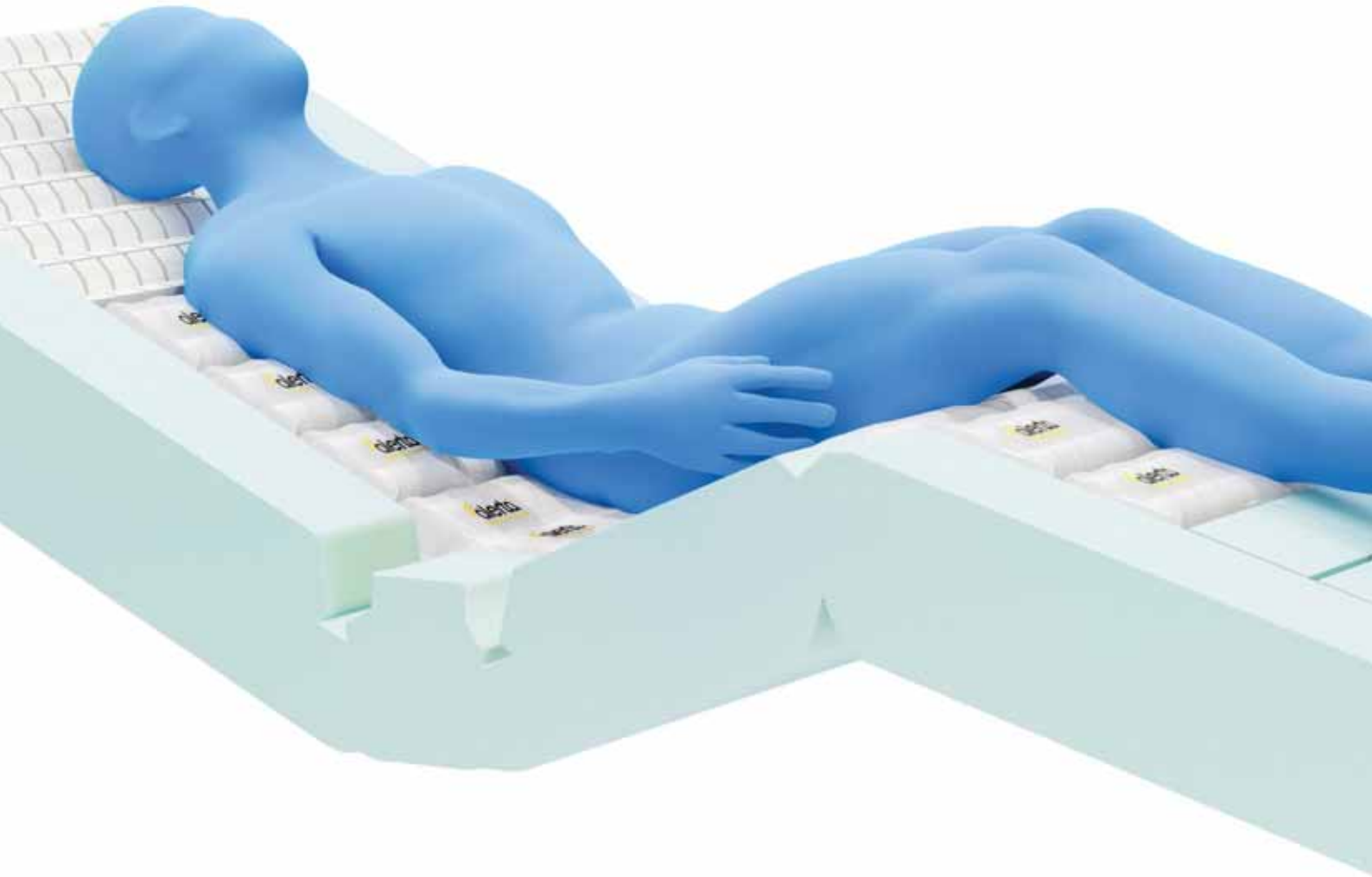



# Alternating Air Systems


# Product Overview

		Overlay pad	
		Overlay system	
		Overlay system	
		Overlay system	
		Replacement system	
		Full Replacement system	
		Full Replacement system	
		Full Replacement system	
		Full Replacement system	
		Cushion system	
		Battery cushion system	





Extension flaps

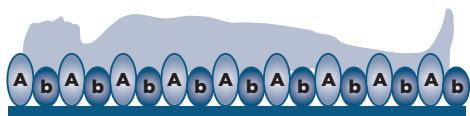

with tuck-in extension flaps

ALT-9001  
Alerta Bubble 2



**Alerta Bubble 2** – a light weight and compact overlay alternating pressure relieving mattress system for effective prevention and treatment of users at low risk of developing a pressure ulcer in hospital, nursing and care home environments. **Alerta Bubble 2** has been designed to give care providers the basic essentials for pressure care and the pump is a basic entry level pump.

With simple to use settings and functionality, the **Alerta Bubble 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The **Alerta Bubble 2** has been manufactured to comply with the most stringent quality and in-use guidelines.


\*If putting in washing machine, ensure the air inlet tubes are watertight sealed.

**Mattress Specification**

- > Mattress Dimensions: 200 x 90 x 7 cm (78.5 x 35.5 x 3")
- > Mattress Weight: 2.5kg (5.5lb)
- > Cell Structure: 130 Bubble cells
- > Material: Medical grade PVC
- > Max Load: 135kg (297.5lb)
- > Machine washable\*

**Pump Specification**

- > Pump Dimensions: 24 x 11 x 8 cm (9.4 x 4.3 x 3.1")
- > Pump Weight: 0.9kg (1.7lb)
- > Cycle Time: 12 minutes
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets

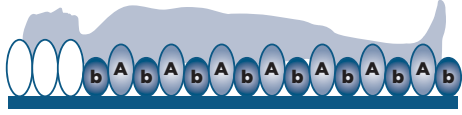

**ALT-9004**  
**Alerta Pearl**  
**Overlay System**



**ALT-9004 - Alerta Pearl** – an overlay alternating pressure relieving mattress system with 5" air cells for effective prevention and treatment of users at medium to high risk of developing a pressure ulcer in hospital, nursing and care home environments. **Alerta Pearl** has been designed to give care providers the essentials for pressure care and the pump is an entry level pump.

With simple to use settings and functionality, the **Alerta Pearl** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The **Alerta Pearl** has been manufactured to comply with the most stringent quality and in-use guidelines.

**Key Features**

- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Micro-perforated cells for air circulation
- > Robust pump connection
- > Transport mode

**Mattress Specification**

- > Mattress Dimensions: 200 x 90 x 12.7 cm (78.5 x 35.5 x 5")
- > Mattress Weight: 5.5kg (12.1lb)
- > CPR: Quick release toggle
- > Max Load: 160kg (353lb)


**Cover Specification**

- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

**Pump Specification**

- > Pump Dimensions: 28 x 14 x 9 cm (11 x 5.5 x 3.5")
- > Pump Weight: 1.4kg (3lb)
- > Cycle Time: 12 minutes
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure and power failure

## KEY FEATURES AND BENEFITS

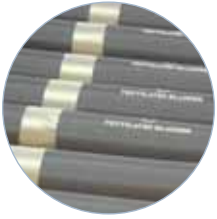

Nylon TPU coated cells




Independently removable cells



Covert sealed air cells



Branded cells

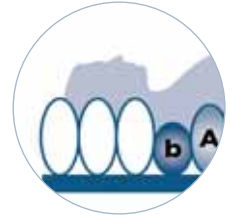

Pillow function



CPR release




Anchoring straps



Zippered storage bag



Multi-stretch, waterproof & vapour permeable PU cover

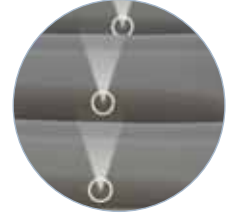

Micro-perforated cells



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode


## CONTROL PANEL

Visual and audible warning of air loss and power failure




Weight setting control knob

Static/alternating modes

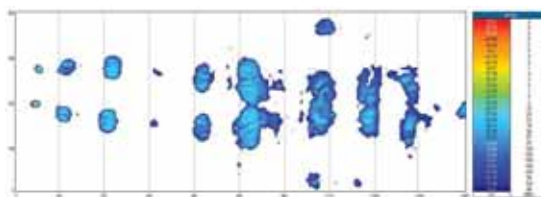

Clinical pressure mapping evaluation

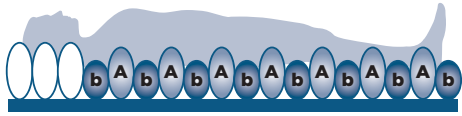

**ALT-9005**  
Alerta Emerald 2  
Overlay System



**ALT-9005 - Alerta Emerald 2** – an overlay alternating pressure relieving mattress system with 5" air cells for effective prevention and treatment of users at high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Emerald 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The wide range of features make the **Alerta Emerald 2** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

**Alerta Emerald 2** incorporates cutting-edge **SensaCare TruFlo Technology®** - constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

**Key Features**

- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Micro-perforated cells for air circulation
- > Robust pump connection
- > Transport mode

**Mattress Specification**

- > Mattress Dimensions: 200 x 90 x 12.7 cm (78.5 x 35.5 x 5")
- > Mattress Weight: 5.5kg (12.1lb)
- > CPR: Quick release toggle
- > Max Load: 160kg (353lb)

**Cover Specification**


- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

**Pump Specification**

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode



## KEY FEATURES AND BENEFITS

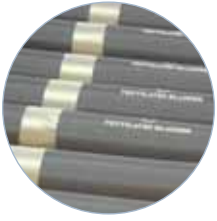

Nylon TPU coated cells




Independently removable cells



Covert sealed air cells




Branded cells

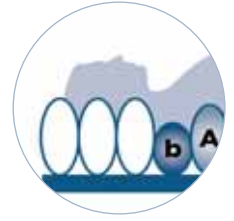

Pillow function




CPR release

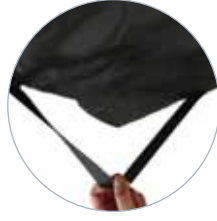


Anchoring straps



Zipped storage bag

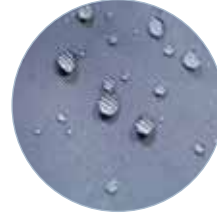

Multi-stretch, waterproof & vapour permeable PU cover

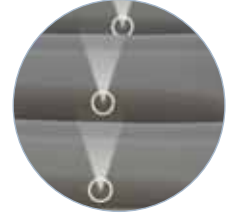

Micro-perforated cells



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode

## CONTROL PANEL

Auto pressure sensing for consistent distribution of air


Visual and audible warning of air loss, power failure and alternate failure

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)




Ultra-Silent™ pump operation to avoid disturbance

Multiple cycle time control for efficiency

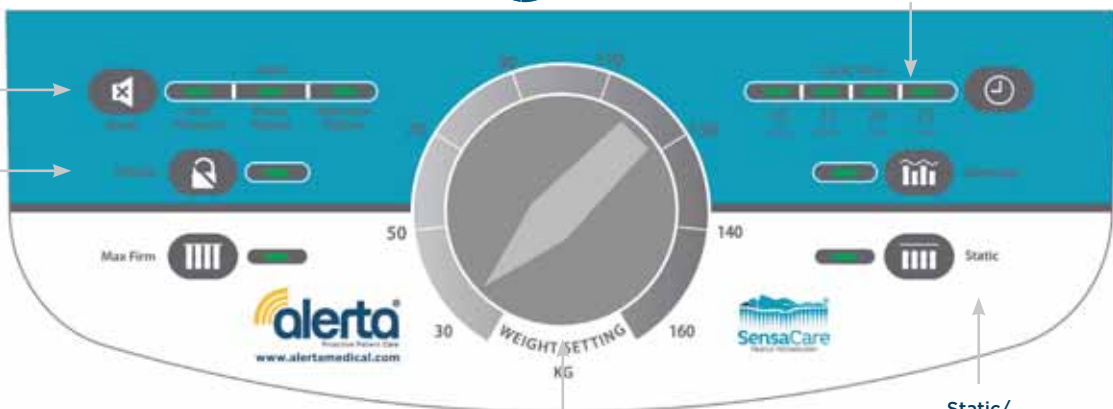


Weight setting control knob

Static/alternating modes

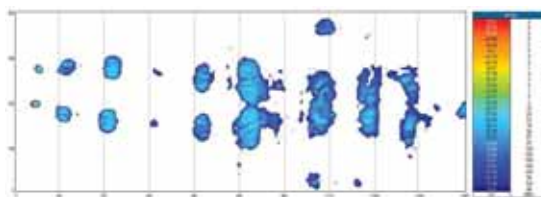

Clinical pressure mapping evaluation

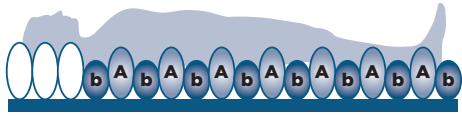

with remote

**ALT-9006**  
Alerta Emerald Auto Overlay System



**ALT-9006 - Alerta Emerald Auto** – a fully **Auto Weight Sensing** overlay alternating pressure relieving mattress system with 5” air cells for effective prevention and treatment of users at high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Emerald Auto** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease. The pump has the latest cutting-edge **Auto Weight Sensing** technology, continuously adjusting interface pressures according to the user’s weight.

The wide range of features make the **Alerta Emerald Auto** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

**Alerta Emerald Auto** incorporates cutting-edge **SensaCare TruFlo Technology®** – constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

**Key Features**

- > Auto weight sensing
- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Robust pump connection
- > Transport mode

**Mattress Specification**

- > Mattress Dimensions: 200 x 90 x 12.7 cm (78.5 x 35.5 x 5")
- > Mattress Weight: 5.5kg (12.1lb)
- > CPR: Quick release toggle
- > Max Load: 160kg (353lb)


**Cover Specification**

- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

**Pump Specification**

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

## KEY FEATURES AND BENEFITS

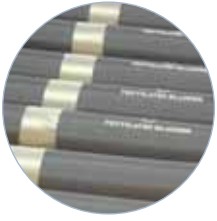

Nylon TPU coated cells




Independently removable cells



Covert sealed air cells




Branded cells

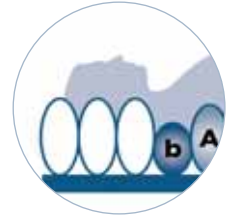

Pillow function



CPR release

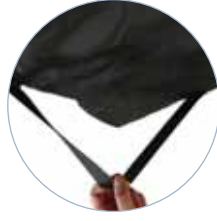

Anchoring straps



Zipped storage bag



Multi-stretch, waterproof & vapour permeable PU cover



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



Auto weight



Remote control

## CONTROL PANEL

Auto pressure sensing for consistent distribution of air



Ultra-Silent™ pump operation to avoid disturbance

Multiple cycle time control for efficiency

Auto weight control

Remote control communication

Visual and audible warning of air loss, power failure and alternate failure

Alternate

Static

Max Firm

Unlock


Static/alternating modes

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)

Power on


Alerta Medical logo and website: www.alertamedical.com

SensaCare logo

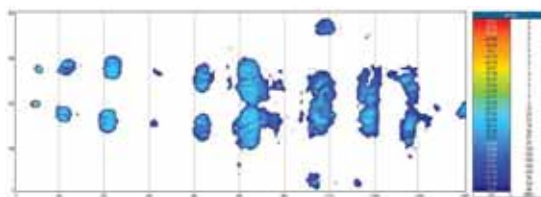

Clinical pressure mapping evaluation

# alerta<sup>®</sup> sapphire 2

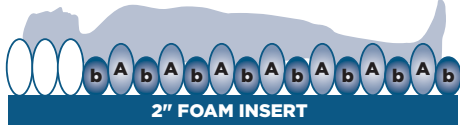

ALT-9007

## Alerta Sapphire 2 Replacement System



**ALT-9007 – Alerta Sapphire 2** – a replacement alternating pressure relieving mattress system with 5" air cells and a self-inflating 2" sealed foam base forming an in-use height of 7" for effective prevention and treatment of users at high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Sapphire 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The wide range of features make the **Alerta Sapphire 2** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

**Alerta Sapphire 2** incorporates cutting-edge **SensaCare TruFlo Technology<sup>®</sup>** - constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

### Key Features

- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Micro-perforated cells for air circulation
- > Robust pump connection
- > Transport mode
- > Individually removable heel cells to create zero pressure under the heel

### Mattress Specification

- > Mattress Dimensions: 200 x 90 x 12.7 cm + 5cm foam base (78.5 x 35.5 x 5" + 2" foam base)
- > Mattress Weight: 5.5kg (12.1lb)
- > Cell Structure: Cell-on-foam
- > CPR: Quick release toggle
- > Max Load: 200kg (441lb)
- > Heel function (3 x heel cells with one way valves)

### Cover Specification


- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

### Pump Specification

- > Operating System: SensaCare TruFlo Technology<sup>®</sup>
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode



## KEY FEATURES AND BENEFITS

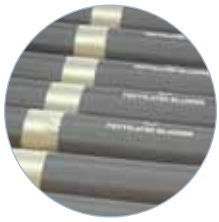

Nylon TPU coated cells




Independently removable cells



Covert sealed air cells



Branded cells

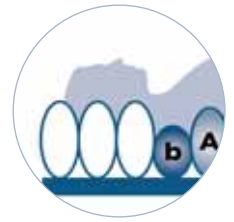

Pillow function




CPR release




Anchoring straps



Zipped storage bag

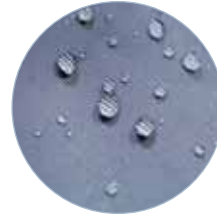

Multi-stretch, waterproof & vapour permeable PU cover

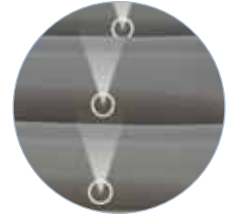

Micro-perforated cells



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



Removable heel cells

## CONTROL PANEL

Auto pressure sensing for consistent distribution of air


Visual and audible warning of air loss, power failure and alternate failure

Control panel lock to prevent unauthorised Or accidental tampering (hold for 2 sec)




Ultra-Silent™ pump operation to avoid disturbance

Multiple cycle time control for efficiency

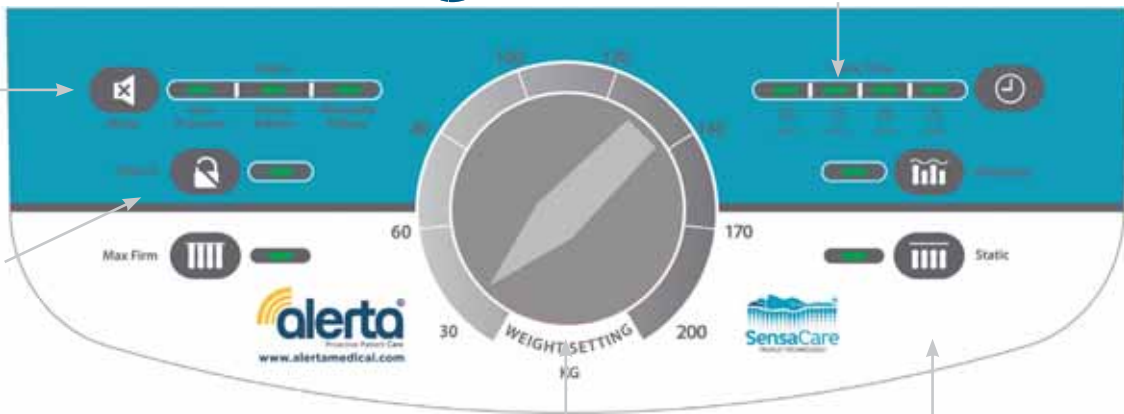


Weight setting control knob

Static/alternating modes

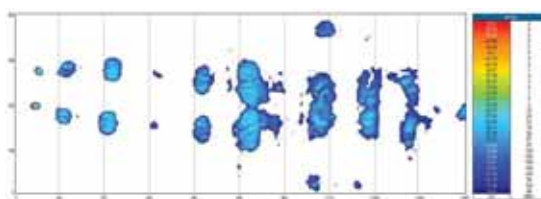

Clinical pressure mapping evaluation

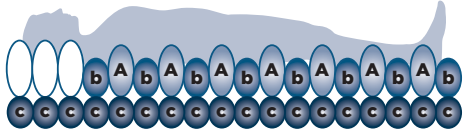

with remote

**ALT-9008**  
**Alerta Ruby 2**  
**Replacement System**

**ALT-9008 - Alerta Ruby 2** - a replacement alternating pressure relieving mattress system with cell-on-cell air cells, a static base section and an in-use height of 8" for effective prevention and treatment of users at very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Ruby 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The wide range of features make the **Alerta Ruby 2** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

**Alerta Ruby 2** incorporates cutting-edge **SensaCare TruFlo Technology®** - constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.



**Key Features**

- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Micro-perforated cells for air circulation
- > Robust pump connection
- > Transport mode
- > Individually removable heel cells to create zero pressure under the heel
- > Seat inflation feature for additional support when sitting upright
- > Digital user weight input

**Mattress Specification**

- > Mattress Dimensions: 200 x 90 x 20 cm (78.5 x 35.5 x 8")
- > Mattress Weight: 10kg (22lb)
- > Cell Structure: Cell-on-cell
- > CPR: Quick release toggle
- > Max Load: 200kg (441lb)
- > Heel function (3 x heel cells with one way valves)


**Cover Specification**

- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

**Pump Specification**

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

## KEY FEATURES AND BENEFITS

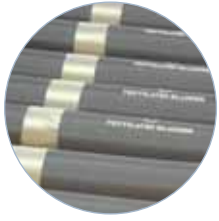

Nylon TPU coated cells



Independently removable cells



Removable heel cells




Branded cells



Pillow function




CPR release

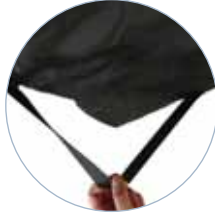


Anchoring straps



Zipped storage bag

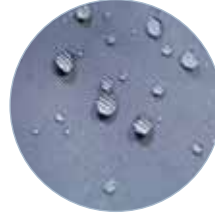

Multi-stretch, waterproof & vapour permeable PU cover

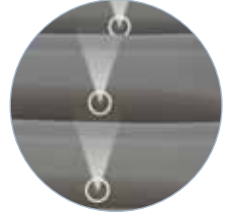

Micro-perforated cells



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



Remote control

## CONTROL PANEL

Digital LED display for easy use

Digital user weight input

Multiple cycle time control for efficiency



Ultra-Silent™ pump operation to avoid disturbance

Visual and audible warning of air loss, power failure and alternate failure

Remote control communication




Static/alternating modes


Seat inflation for additional surface support to facilitate the sitting position

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)

Power on

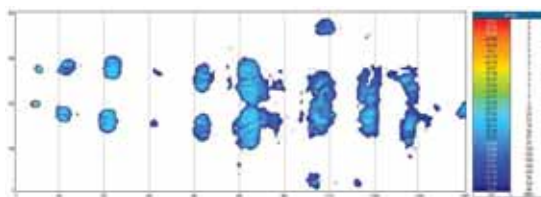

Clinical pressure mapping evaluation

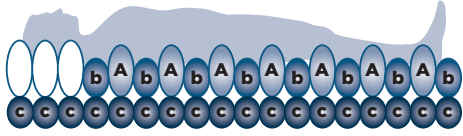

with remote

**ALT-9009**  
Alerta Ruby Auto Replacement System



**ALT-9009 - Alerta Ruby Auto** – a fully **Auto Weight Sensing** replacement alternating pressure relieving mattress system with cell-on-cell air cells, a static base section and an in-use height of 8” for effective prevention and treatment of users at very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Ruby Auto** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease. The pump has the latest cutting-edge **Auto Weight Sensing** technology, continuously adjusting interface pressures according to the user’s weight.

The wide range of features make the **Alerta Ruby Auto** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

**Alerta Ruby Auto** incorporates cutting-edge **SensaCare TruFlo Technology®** – constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

**Key Features**

- > Auto weight sensing
- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Robust pump connection
- > Transport mode
- > Individually removable heel cells to create zero pressure under the heel
- > Seat inflation feature for additional support when sitting upright

**Mattress Specification**

- > Mattress Dimensions: 200 x 90 x 20 cm (78.5 x 35.5 x 8")
- > Mattress Weight: 10kg (22lb)
- > Cell Structure: Cell-on-cell
- > CPR: Quick release toggle
- > Max Load: 200kg (441lb)
- > Heel function (3 x heel cells with one way valves)

**Cover Specification**


- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

**Pump Specification**

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode



## KEY FEATURES AND BENEFITS

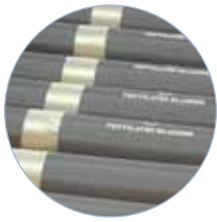

Nylon TPU coated cells



Independently removable cells



Removable heel cells



Branded cells



Pillow function




CPR release



Anchoring straps



Zipped storage bag

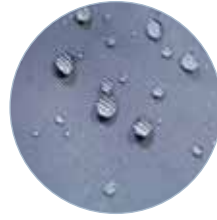

Multi-stretch, waterproof & vapour permeable PU cover



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



Auto weight



Remote control

## CONTROL PANEL

Digital LED display for easy use



Ultra-Silent™ pump operation to avoid disturbance

Multiple cycle time control for efficiency

Auto weight control

Remote control communication

Visual and audible warning of air loss, power failure and alternate failure

Alerta [www.alertamedical.com](http://www.alertamedical.com)

Alternate

Static

Seat Inflate

Max Firm

Unlock


SensaCare

Power on


Static/alternating modes

Seat inflation for additional surface support to facilitate the sitting position

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)

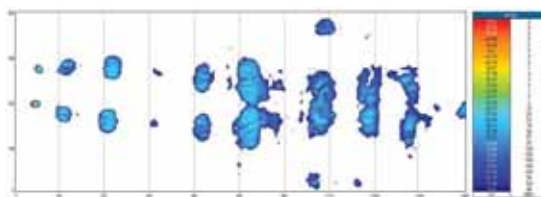

Clinical pressure mapping evaluation

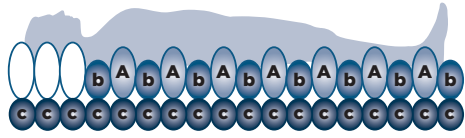

with remote

**ALT-9300**  
**Alerta Bariatric 2**  
**Replacement System**



**ALT-9300 – Alerta Bariatric 2** – a replacement alternating pressure relieving mattress system with cell-on-cell air cells, a static base section and an in-use height of 10” for effective prevention and treatment of bariatric users up to very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Bariatric 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The wide range of features make the **Alerta Bariatric 2** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

**Alerta Bariatric 2** incorporates cutting-edge **SensaCare TruFlo Technology®** - constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

**Key Features**

- > 350kg weight capacity
- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Robust pump connection
- > Transport mode
- > Individually removable heel cells to create zero pressure under the heel
- > Seat inflation feature for additional support when sitting upright
- > Digital user weight input

**Mattress Specification**

- > Mattress Dimensions: 200 x 120 x 25 cm (78.5 x 47 x 10”)
- > Mattress Weight: 12.5kg (27.6lb)
- > Cell Structure: Cell-on-cell
- > CPR: Quick release toggle
- > Max Load: 350kg (772lb)
- > Heel function (3 x heel cells with one way valves)


**Cover Specification**

- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

**Pump Specification**

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9”)
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 12-13 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 9W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

## KEY FEATURES AND BENEFITS

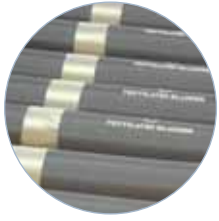

Nylon TPU coated cells



Independently removable cells



Covert sealed air cells




Removable heel cells



Pillow function



CPR release

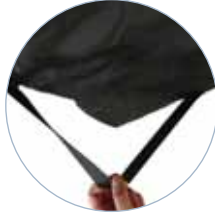

Anchoring straps



Multi-stretch, waterproof & vapour permeable PU cover



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



350KG Capacity



Remote control

## CONTROL PANEL

Digital LED display for easy use

Digital user weight input

Multiple cycle time control for efficiency



Ultra-Silent™ pump operation to avoid disturbance

Visual and audible warning of air loss, power failure and alternate failure

Remote control communication




Static/alternating modes

Seat inflation for additional surface support to facilitate the sitting position

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)


Power on

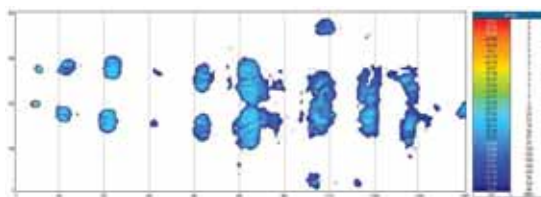

Clinical pressure mapping evaluation



**BARIATRIC  
120CM VERSION**  
Available in 90cm  
and 120cm

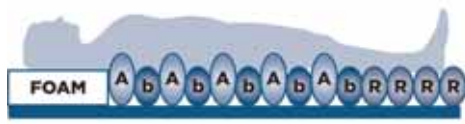

**ALT-HYB2**  
Alerta Sensaflo Hybrid  
mattress only

**ALT-HYB2/4**  
Bariatric 4ft wide Alerta  
Sensaflo Hybrid mattress only

**ALT-HYB2/05**  
Alerta Sensaflo Hybrid  
mattress with pump

**ALT-HYB2/05/4**  
Bariatric 4ft wide Alerta Sensaflo  
Hybrid mattress with pump



The **Alerta Sensaflo Hybrid** is the world's first 4 in 1 hybrid mattress with foam-filled air cells, SensaGel adaptive foot cells and 150mm (6") deep air cells for effective prevention and treatment of users up to very high risk of developing a pressure ulcer.

When used without a pump the **Alerta Sensaflo Hybrid** performs as a high specification foam mattress and when used with a pump it performs the same as an alternating air system. There is no requirement to remove the user from bed to change from a foam system to an alternating air system.

The **Alerta Sensaflo Hybrid** has a built-in fire evacuation system with straps and handles, allowing the mattress to be quickly moved with the user in position during emergency situations.

There are four adaptable gel-foam heel zone cells that can be individually removed and repositioned to create zero pressure under the heel. By inflating the mattress and sealing off with the transport cap you can create a balanced-air self-adjusting semi-active mattress.

The **Alerta Sensaflo Hybrid** incorporates cutting-edge **SensaCare Truflo Technology®**; with constant pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure, and ultra-silent pump operation to promote an even ambience for natural sleep patterns and stable blood circulation.






WATCH PRODUCT DEMO VIDEO

Ergonomically castellated pillow function

In the centre of the mattress are 10 foam filled air cells which the user interfaces with directly

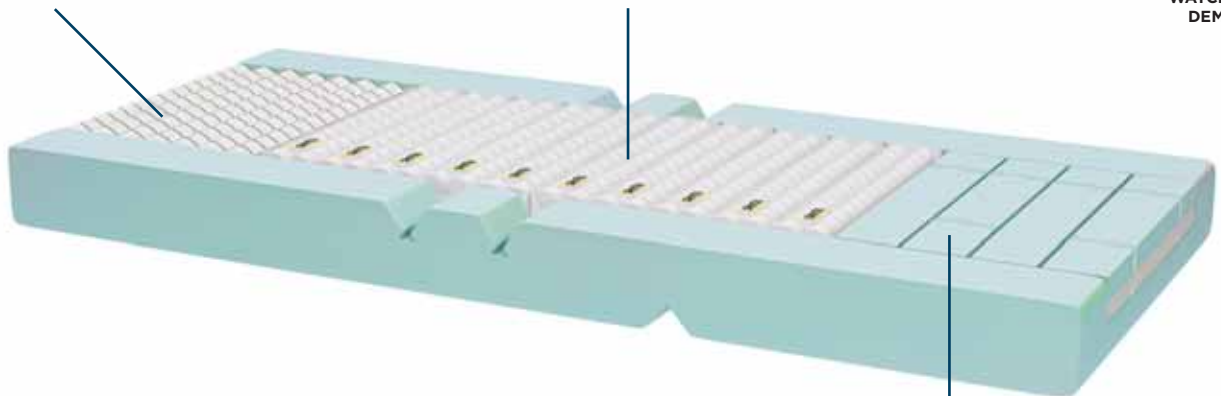

Four adaptable gel-foam heel zone cells can be individually removed and repositioned to create zero pressure under the heel



Built-in fire evacuation system with straps and handles, allow the mattress to be quickly moved with the user in position during emergency situations

### Key Features

- > Pillow function
- > CPR release
- > Robust pump connection
- > Transport mode
- > Profiling side sections
- > Adaptable gel-foam heel zone
- > Built-in evacuation system

### Mattress Specification

- > For users up to very high risk of developing a pressure ulcer
- > CPR quick release toggle
- > Max load: 250kg
- > Min load: 30kg
- > Mattress depth: 150mm (6")

### Cover Specification

- > Multi-stretch and vapour permeable PU cover
- > Water resistant and machine washable cover
- > Antimicrobial cover with white underside
- > 360 degree zipper
- > BS7177: Crib 5
- > Mains cable management system

### Pump Specification

- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control Panel lock feature to avoid accidental change of settings
- > Fitted with hanging brackets
- > Max firm mode

# alerta<sup>®</sup> sensafllo-hybrid



## ALT-HYB-C-SYS

### Alerta Sensafllo Hybrid cushion with pump



**ALT-HYB-C-SYS - Alerta Sensafllo Hybrid cushion system** – an alternating foam-filled air cell pressure relieving cushion system for effective prevention and treatment of users up to very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Sensafllo Hybrid** cushion enables the care provider to quickly set up the cushion system and have it operating at optimal pressure for a specific user with ease. There is no requirement to remove the user from a chair to change from a foam cushion to an alternating air system.

**Alerta Sensafllo Hybrid** cushion can be used without a pump and boasts a unique SensaGel gel-infused memory foam topper for unrivalled pressure redistribution and improved circulation. When used with a pump it performs the same as an alternating air system. By inflating the cushion and sealing off with the transport cap you can create a balanced-air self-adjusting semi-active cushion.

The **Alerta Sensafllo Hybrid** incorporates cutting-edge **SensaCare Truflo Technology<sup>®</sup>**; with constant pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure, and ultra-silent pump operation to promote an even ambience for natural comfort and stable blood circulation.

### Cushion Specification


- > Cushion Dimensions:  
45x45x12cm inflated including SensaGel insert  
45x45x10cm uninflated  
45x45x8cm uninflated & SensaGel insert removed
- > Cushion Weight: 1.5kg (3.3lb)
- > Cell Structure: Foam filled cell
- > Max Load: 160kg (353lb)

### Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5

### Pump Specification

- > Pump Dimensions:  
27 x 15.5 x 10 cm (10.5 x 6.1 x 3.9")
- > Pump Weight: 1.6kg (3.5lb)
- > Cycle Time: 10 minutes
- > Power Input:  
AC220-240V, 50-60Hz
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure


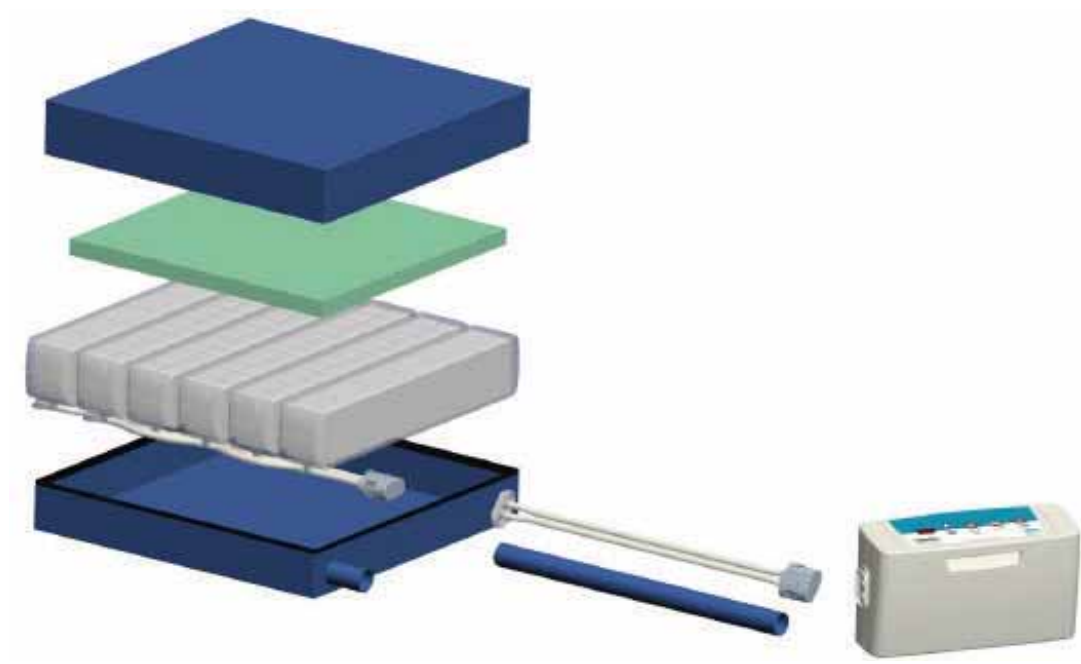


**CONTROL PANEL**  
Digital LED display for easy use

Digital user weight input



Ultra-Silent™ pump operation to avoid disturbance

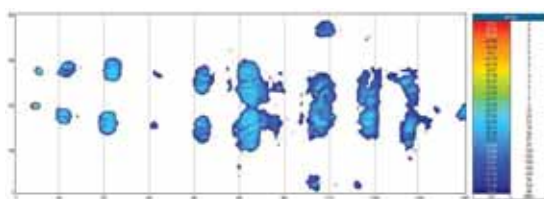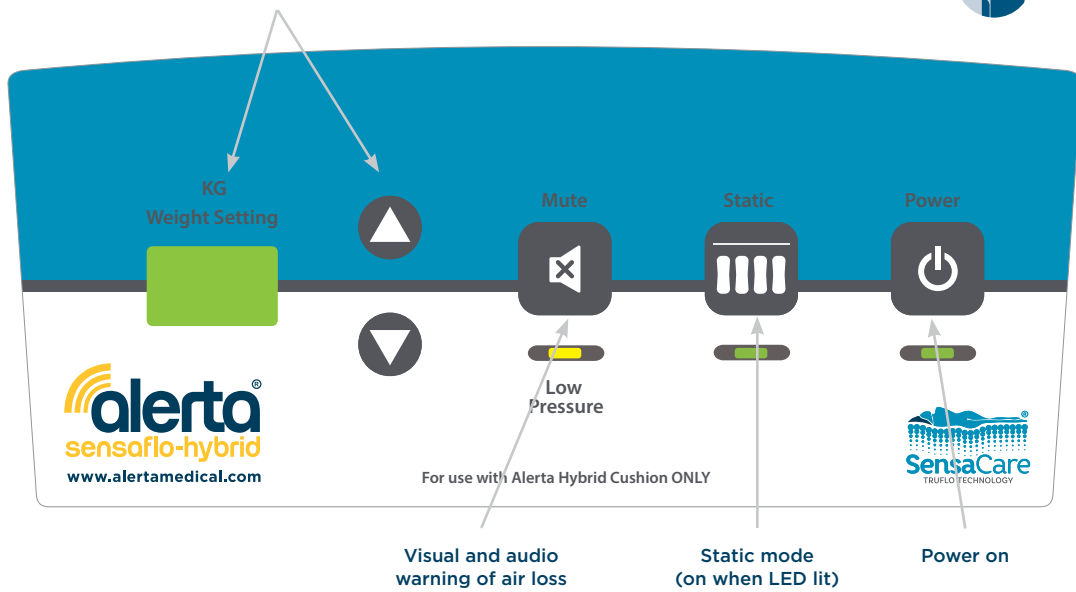

Clinical pressure mapping evaluation




**ALT-213/02**  
Alerta Partner  
Cushion System

**ALT-213/02 – Alerta Partner** – an alternating pressure relieving cushion system with a sealed 2" foam base for effective prevention and treatment of users up to high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Partner** enables the care provider to quickly set up the cushion system and have it operating at optimal pressure for a specific user with ease.

The **Alerta Partner** is a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

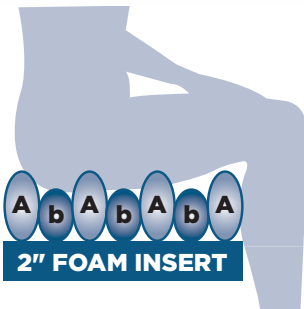

**Cushion Specification**

- > Cushion Dimensions:  
45 x 50 x 10 cm  
(17.5 x 19.5 x 4")
- > Cushion Weight: 1.4kg (3.1lb)
- > Cell Structure: Cell-on-foam
- > Max Load: 150kg (331lb)


**Pump Specification**

- > Pump Dimensions:  
27 x 13.5 x 9 cm (10.5 x 5 x 3.5")
- > Pump Weight: 1.6kg (3.5lb)
- > Cycle Time: 12 minutes
- > Power Input:  
AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure and power failure

**Cover Specification**

- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5

**CONTROL PANEL**

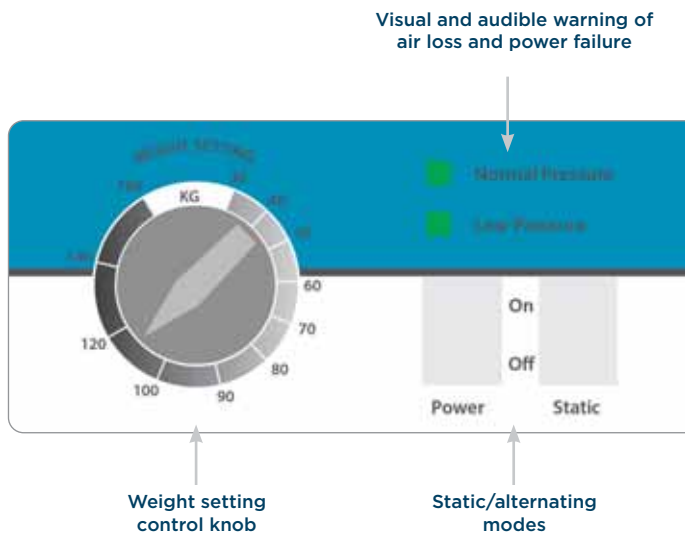

**KEY FEATURES AND BENEFITS**




Nylon TPU coated cells



Independently removable cells



Multi-stretch, waterproof & vapour permeable PU cover

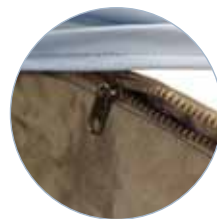

360 Degree zipper



Robust pump connection



Transport mode






**ALT-215/04**  
Alerta Mobile Battery  
Cushion System



**ALT-215/04 – Alerta Mobile** – a battery powered alternating pressure relieving cushion system for effective prevention and treatment of users up to high risk of developing a pressure ulcer in hospital, nursing and care home environments. The **Alerta Mobile** control unit is equipped with a long lasting Lithium Ion battery. In static mode the battery gives a remote operation time of more than 24 hours and in alternating mode the system can last at least one day.

With simple to use settings and functionality, the **Alerta Mobile** enables the care provider to quickly set up the cushion system and have it operating at optimal pressure for a specific user with ease.

The **Alerta Mobile** is a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

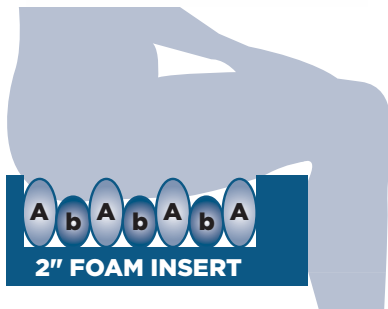

**Cushion Specification**

- > Cushion Dimensions: 45 x 50 x 10 cm (17.5 x 19.5 x 4")
- > Cushion Weight: 2.4kg (5.3lb)
- > Cell Structure: Cell-on-foam
- > Max Load: 130kg (286.6lb)


**Pump Specification**

- > Pump Dimensions: 21 x 10 x 6 cm (8.2 x 4 x 2.5")
- > Pump Weight: 0.7Kg (1.54lb)
- > Cycle Time: 10, 15, 20 minutes
- > Power Input: AC220-240V, 50-60Hz
- > Battery: 12V Lithium 2200mAh
- > Enclosed pump case fitted with hanging straps
- > Alternating and static modes
- > Visual and audible alarm for low pressure and low battery power


**Cover Specification**

- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > Enclosed zipper
- > White underside of cover
- > BS7177: Crib 5

**CONTROL PANEL**

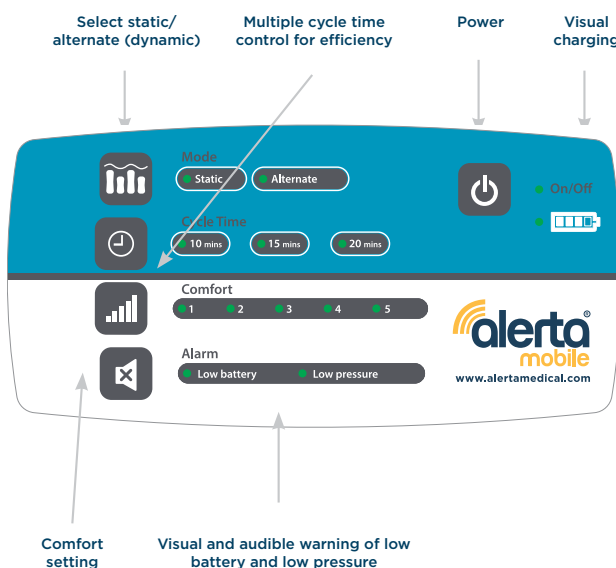

**KEY FEATURES AND BENEFITS**

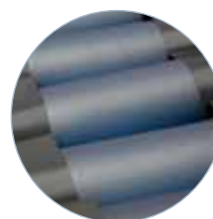

TPU coated cells



Independently removable cells



Multi-stretch, waterproof & vapour permeable PU cover













Antimicrobial, machine washable



360 Degree zipper

# Product Selector

										
User Risk Level	High/Very High	Very High	Very High	Very High	High	High	High	Med/High	High	High
Max User Weight (kg)	250	350	200	200	200	160	160	160	150	130
Cell structure	Foam and air/foam	Cell-on-cell	Cell-on-cell	Cell-on-cell	Cell-on-cell -on-foam	Cell-on-cell	Cell-on-cell	Cell-on-cell	Cell-on-cell -on-foam	Cell-on-cell -on-foam
Overlay system						✓	✓	✓		
Replacement system	✓	✓	✓	✓	✓				✓	✓
Auto weight sensing			✓			✓				
Built-in evacuation system	✓									
Nylon TPU coated cells		✓	✓	✓	✓	✓	✓	✓	✓	✓
Independently removable cells	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
Heel zone	✓		✓	✓	✓					
Covert sealed air cells		✓	✓	✓	✓	✓	✓	✓	✓	✓
Branded cells	✓	✓	✓	✓	✓	✓	✓	✓		
Pillow function	✓	✓	✓	✓	✓	✓	✓	✓		
CPR	✓	✓	✓	✓	✓	✓	✓	✓		
Anchoring straps		✓	✓	✓	✓	✓	✓	✓		
Zippered storage bag		✓	✓	✓	✓	✓	✓	✓		
Multi-stretch PU cover	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
Waterproof, vapour permeable	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
Antimicrobial, machine washable	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
360 degree zipper	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
Robust pump connection	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
Transport mode	✓	✓	✓	✓	✓	✓	✓	✓		
Remote control		✓	✓	✓		✓				
Clinical pressure mapping evaluation	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓

# Bespoke Branding

*"Alerta Medical offer a fast and cost effective own branding solution for alternating mattress systems globally with no MOQ"*

**Alerta Medical offer a bespoke branding service at no extra cost and with no minimum order quantity - truly enhancing your brand and increasing your profitability.**

As a distributor of Alerta Medical's alternating air mattress systems you can have your company logo and contact details prominently displayed on these products - helping to enhance your brand and set yourself apart from the competition.

**Be proactive - contact Alerta Medical now to take advantage of this unique bespoke branding service!**



## MANUFACTURER PRICES

As the first step in the supply chain, we offer an attractive pricing structure

**- increasing your profitability as a distributor**




## OWN BRANDED PRODUCTS

At no extra cost and no minimum order.



## CUSTOMISED STOCKHOLDING

Guaranteed bespoke branded stock based on pre-determined requirements

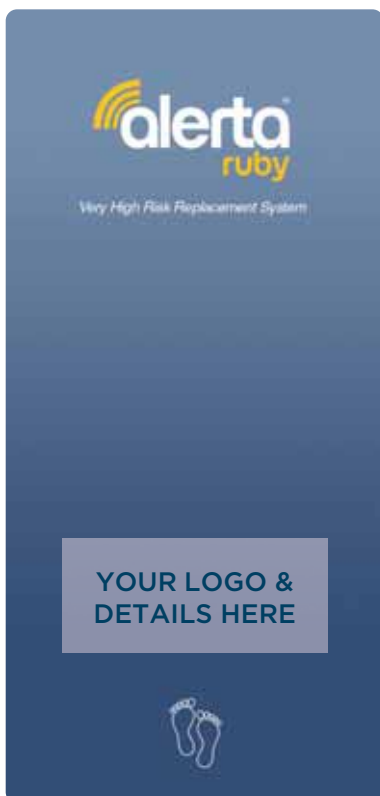

**The bespoke branding service is available on the following products:**

- > Alerta Pearl
- > Alerta Emerald 2
- > Alerta Emerald Auto
- > Alerta Sapphire 2
- > Alerta Ruby 2
- > Alerta Ruby Auto
- > Alerta Bariatric 2

# Suction Machines

## AS1003A

### 3A PROFESSIONAL ASPIRATOR SINGLE PUMP, 1000ML BOTTLE



- > 16 litre/minute flow rate
- > Adjustable vacuum level
- > Machine only weighs 2.3kg
- > Carry handle

#### What's included?

- > Machine
- > 1000ml jar
- > Suction tubing
- > 1 suction catheter
- > Spare anti-bacterial filter
- > Mains power lead
- > Connector/flow regulator
- > Disposable jar liner

## AS903A

### 3A PROFESSIONAL ASPIRATOR DOUBLE PUMP, 1000ML BOTTLE



- > 22 litre/minute flow rate
- > Adjustable vacuum level
- > Machine only weighs 3.2kg
- > Carry handle

#### What's included?

- > Machine
- > 1000ml jar
- > Suction tubing
- > 1 suction catheter
- > Spare anti-bacterial filter
- > Mains power lead
- > Connector/flow regulator
- > Disposable jar liner

## AS457A

### 3A MINI RECHARGEABLE PORTABLE ASPIRATOR



- > 30 litre/minute flow rate
- > Adjustable vacuum level
- > Machine only weighs 4.6kg
- > Carry handle
- > Rechargeable for independent use
- > Can also be mains operated
- > 12V connection facility

#### What's included?

- > Machine
- > 1000ml jar
- > Suction tubing
- > 1 suction catheter
- > Spare anti-bacterial filter
- > Mains power lead
- > Connector/flow regulator
- > Disposable jar liner
- > Shoulder bag/carry case



# Nebulisers

## RJ-2 NEBULISER



- > 0.26 ml/min nebulisation rate
- > Mains operated

### What's included?

- > Nebuliser
- > Power cable
- > Tubing
- > Face mask
- > Mouthpiece
- > Filter

## AP508A HAPPYNEB III NEBULISER



- > 14 litre/minute flow rate
- > .45ml/minute nebulisation rate
- > Machine weighs 2kg
- > Mains operated

### What's included?

- > Nebuliser
- > Power cable
- > Tubing
- > 2 face masks
- > Chamber
- > Mouth piece
- > Nose piece
- > Filter

# Accessories

## 3A1385 ANTIBACTERIAL FILTER FOR3A PROFESSIONAL ASPIRATORS



## 3A561 SILICON SUCTION TUBING 1.3M LONG



## 3A560 MANUAL FLOW REGULATOR & CONNECTOR



## 3AC285 REPLACEMENT VASE 1000ML WITH CONVERSION KIT



## 3A1687 DISPOSABLE SUCTION LINER



## 3AC170 FACE MASK KIT FOR NEBULISER

### What's included?

- > Face mask
- > Air tube
- > Chamber



## 3AC174 POREX FILTER FOR AP103A NEBULISER



# Diagnostics

## ALT-DE100

### DIGITAL NON-CONTACT FOREHEAD & EAR THERMOMETER

An easy-to-use functional thermometer for recording body temperature quickly and accurately. It holds a memory of the last 30 readings and automatically shuts down after 15 seconds of no use.



## ALT-DE102

### DIGITAL BODY THERMOMETER

An easy-to-use functional thermometer for recording body temperature quickly and accurately. It holds a memory of the last reading and automatically shuts down after 10 minutes of no use.



## ALT-DE200

### UPPER ARM DIGITAL BLOOD PRESSURE MONITOR

An easy-to-use functional device for recording blood pressure quickly and accurately. It measures systolic pressure, diastolic pressure and pulse. With automatic inflation and deflation, this device is comfortable to wear and provides a clinically precise reading in seconds.

It has a memory capacity for 2 sets of 90 readings, provides an average of the latest 3 readings and automatically shuts down after no use.

A cuff is provided and fits upper arm circumferences of 22-32cm.



## ALT-DE300

### FINGERTIP PULSE OXIMETER

An easy-to-use functional oximeter for recording blood oxygen saturation and pulse rate quickly and accurately by simply placing on the finger. This device is comfortable to use and provides a clinically precise reading in seconds. It automatically turns off after no use and has a low battery indicator.





